

# INTEGRA



**Linear Stand**  
**Charging Stations for Linear Stand**

**Operating instructions**



# Declaration of conformity

**INTEGRA Biosciences AG – 7205 Zizers, Switzerland**

declares on its own responsibility that the devices

Description	Models
<b>Linear Stand</b>	<b>3215</b>
<b>Charging Station for Linear Stand</b>	<b>3216</b>
<b>Charging/Communication Station for Linear Stand</b>	<b>3217</b>

comply with:

<b>EU Directives</b> (DoW: Date of Withdrawal)	Before DoW	DoW	After DoW
Low Voltage Equipment	<b>2006/95/EC</b>	20.04.2016	<b>2014/35/EU</b>
Electromagnetic Compatibility	<b>2004/108/EC</b>	20.04.2016	<b>2014/30/EU</b>
Restriction of Hazardous Substances	<b>2011/65/EU</b>		
Waste Electrical and Electronic Equipment	<b>2012/19/EU</b>		

## EU Regulation

Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH) **1907/2006**

## Standards for EU

Safety requirements for electrical equipment for measurement, control and laboratory use - General requirements.	<b>EN 61010-1: 2010</b>
Electrical equipment for measurement, control and laboratory use - EMC requirements.	<b>EN 61326-1: 2013</b>

## Standards for Canada and USA

Safety requirements for electrical equipment for measurement, control and laboratory use - General requirements.	<b>CAN/CSA-C22.2 No. 61010-1</b>
Safety requirements for electrical equipment for measurement, control and laboratory use - General requirements.	<b>UL 61010-1</b>

Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.	<b>Part 15 of the FCC Rules Class A</b>
---	---

Zizers, November 4, 2016

Elmar Morscher  
CEO

Thomas Neher  
Quality Manager

## Table of Contents

1	Introduction .....	3
2	Installation .....	3
3	Disposal .....	5
4	Accessories .....	5

## 1 Introduction

### 1.1 Intended use

The Linear Stand is used to hang up your INTEGRA Biosciences pipettes.

The Charging Station and the Charging/Communication Station, here Charging Stations for short, are mounted on the Linear Stand (#3215) to charge VIAFLO/VOYAGER Electronic Pipettes from INTEGRA Biosciences. Using the Charging/Communication Station, the pipettes can be connected to VIALINK via a PC. VIALINK allows the user to upload/download custom programs, images, firmware updates and service history to and from VIAFLO and VOYAGER pipettes.

The Charging Stations have been designed for use in a laboratory. They shall be used in a dry and dust-free location with a temperature of 5–40 °C and a maximal (non-condensing) relative humidity of 80 %.


### 1.2 Safety notes

- 1) Do not use the Charging Stations in an atmosphere with danger of explosion.
- 2) Use an original INTEGRA Biosciences Mains Adapter (#3216) only.
- 3) Prolonged exposure of Charging Stations to UV-light can cause discoloration and/or yellowing of the plastic housing. This will not affect the performance of the device in any way.

Please visit our website [www.integra-biosciences.com](http://www.integra-biosciences.com) on a regular basis for up to date information regarding REACH classified chemicals contained in our products.


## 2 Installation

### 2.1 Assembling the Linear Stand

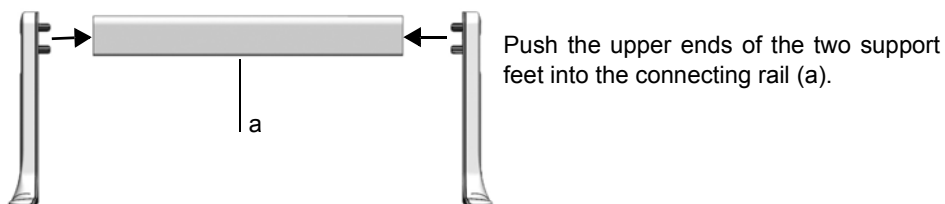


## 2.2 Mounting the Charging Stations

**Note:** The Charging Stations (#3217, #3218), Linear Stand (#3215) and Mains Adapter (#3216) must be ordered separately.

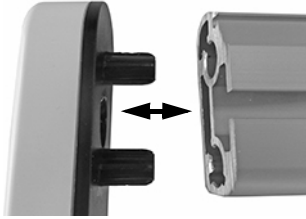


- 1) Take the Linear Stand (#3215) and pull one of the two support feet off the connecting rail.

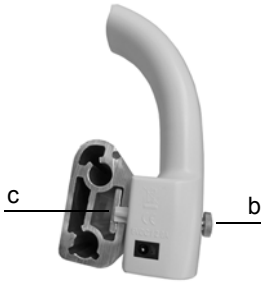


- 2) Loosen the screw (b) on the Charging Station and slide the rectangular washer (c) of the Charging Station on the connecting rail.
- 3) Move the Charging Station to the desired position and tighten the screw. Up to 4 Charging Stations can be mounted to each linear stand.



- 4) Insert the plug (d) of mains adapter into the bottom of the removed support foot and guide it through the foot.

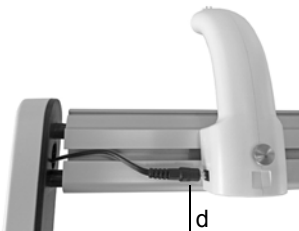

- 5) Pass the cable through the horizontal slot of the connecting rail.
- 6) Connect the plug (d) with the Charging Station.

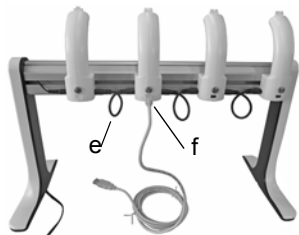

- 7) Push the support foot back on the Linear Stand.
- 8) Connect the Charging Stations serially using the connection cable (e).
- 9) Insert the plug of the USB cable (f) into the socket of the Charging/Communication Station, if applicable.

## 3 Disposal



The Linear Stand can be disposed of with municipal waste. The Charging Stations and the Mains Adapter must not be disposed of with unsorted municipal waste. Do not dispose them in a fire.

Dispose of the Charging Stations and the Mains Adapter in accordance with the laws and regulations in your area governing disposal of devices.

## 4 Accessories

Accessories		Part No.
Linear Stand	holds up to 4 Charging Stations	3215
Mains Adapter	for use with Linear Stand (100–240 VAC, 50/60 Hz)	3216
Charging Station for Linear Stand	for INTEGRA Electronic Pipettes (incl. connection cable)	3217
Charging/Communication Station for Linear Stand	for INTEGRA Electronic Pipettes (incl. connection cable and USB cable)	3218

## Imprint

© 2017 INTEGRA Biosciences AG

### Manufacturer

#### **INTEGRA Biosciences AG**

CH-7205 Zizers, Switzerland

T +41 81 286 95 30

F +41 81 286 95 33

[info@integra-biosciences.com](mailto:info@integra-biosciences.com)

[www.integra-biosciences.com](http://www.integra-biosciences.com)

#### **INTEGRA Biosciences Corp.**

Hudson, NH 03051, USA

T +1 603 578 5800

F +1 603 577 5529