

VIALAB Operating instructions

Document Change History

Vers.	Status:	Name:	Date:	Change:
V00		ULe, SSt	07-12-2017	Document issued
V01	Released	SSt	15-03-2018	Document released
V02	Released	SSt	10-10-2018	download information added, troubleshooting "Program Transfer To Pipette" added
V03	Released	ULe, SSt	09-01-2019	New Software V2.0
V04	Released	SSt	13-02-2019	System Requirements for VIALAB 2.0 changed
V05	Released	SSt	04-03-2020	New Software V2.1 – new features
V06	Released	MFr	21-04-2020	Update Pipette Communication Module
V07	Released	JBa2	13-08-2020	Typo removed

Table of Contents

1	Introduction	4
1.1	Purpose	4
1.2	Scope	4
2	System requirements	4
2.1	Personal computer / Laptop	4
2.2	Compatibility	4
2.3	Programming Stand for electronic pipettes	4
2.4	Pipette Communication Module for VIAFLO/VOYAGER electronic pipettes	4
3	Installation of VIALAB	5
3.1	Preparation	5
3.2	First time installation	5
3.3	Running VIALAB without installation	5
4	Settings	6
4.1	VIALAB Settings and Information	6
4.2	Update VIALAB	7
5	Main Menu	8
5.1	Information and saving options	8
5.1.1	Program Settings	8
6	Programming	10
6.1	Overview	10
6.2	Material	12
6.3	Method	14
6.4	Simulation	21
6.5	Transfer	23
6.5.1	Connect pipette	23
6.5.2	Transfer	23
7	Labware Library	24
7.1	Reservoirs	24
7.2	Plates	26
7.3	Tube Rack	28
7.4	COMBI Base	29
7.5	COMBI System	30
7.5.1	Use of a Slanted Plate Holder (PN 4510) in VIALAB	30

1 Introduction

1.1 Purpose

VIALAB enables creating programs and transfers programs to INTEGRA VIAFLO/VOYAGER electronic pipettes. VIALAB programs for use with ASSIST PLUS only.

Programming in VIALAB can be done without a connected pipette.

1.2 Scope

These operating instructions describe the installation and operation of VIALAB V 2.1 or higher.

2 System requirements

2.1 Personal computer / Laptop

CPU: INTEL or AMD (not older than 5 years)
 USB port: 2.0 / 3.0
 Bluetooth: optional
 Memory: > 4 GB
 Internet access: recommended during installation, required for updates
 Operating system: MS Windows 8.1, MS Windows 10.x
 Framework: Microsoft .NET 4.7.2 (or later)

2.2 Compatibility

Device:	Serial No.:	Pipette Firmware connectable to		
		VIALAB V 1.0.x	VIALAB V 2.0.x	VIALAB V 2.1.x
VIAFLO electronic multichannel pipette	7000000 and higher	4.12 and higher	4.20 and higher	4.21 and higher
VOYAGER electronic tip spacing pipette	7000000 and higher	4.12 and higher	4.20 and higher	4.21 and higher

2.3 Programming Stand for electronic pipettes

The Programming Stand (PN 4211) needs to be connected to the USB port of your PC

It can be used to charge the pipette battery and enables the communication between the PC and the pipette.

2.4 Pipette Communication Module for VIAFLO/VOYAGER electronic pipettes

For a wireless communication the Pipette Communication Module (PN 4221) can be used instead of the USB Programming Stand.

3 Installation of VIALAB

3.1 Preparation

Download VIALAB from www.integra-biosciences.com/download-vialab. Register, download the ZIP file and unzip the file.

3.2 First time installation

Log in to your computer with administrator log in.

- Run the setup routine and choose a directory where VIALAB should be installed (e.g. C:\Programs\Integra Biosciences AG\VIALAB). If you are not logged in as administrator, you may be prompted for the administrator login.
- VIALAB can also be installed on a network drive. Choose your network drive during installation.

3.3 Running VIALAB without installation


VIALAB can be run from a USB drive or any other program folder. Simply copy the installed directory (e.g. C:\Programs\Integra Biosciences AG\VIALAB) to your specific drive.

4 Settings

4.1 VIALAB Settings and Information



Click on the Settings button in the top right corner. Use the VIALAB Settings to define your preferences and get additional information.

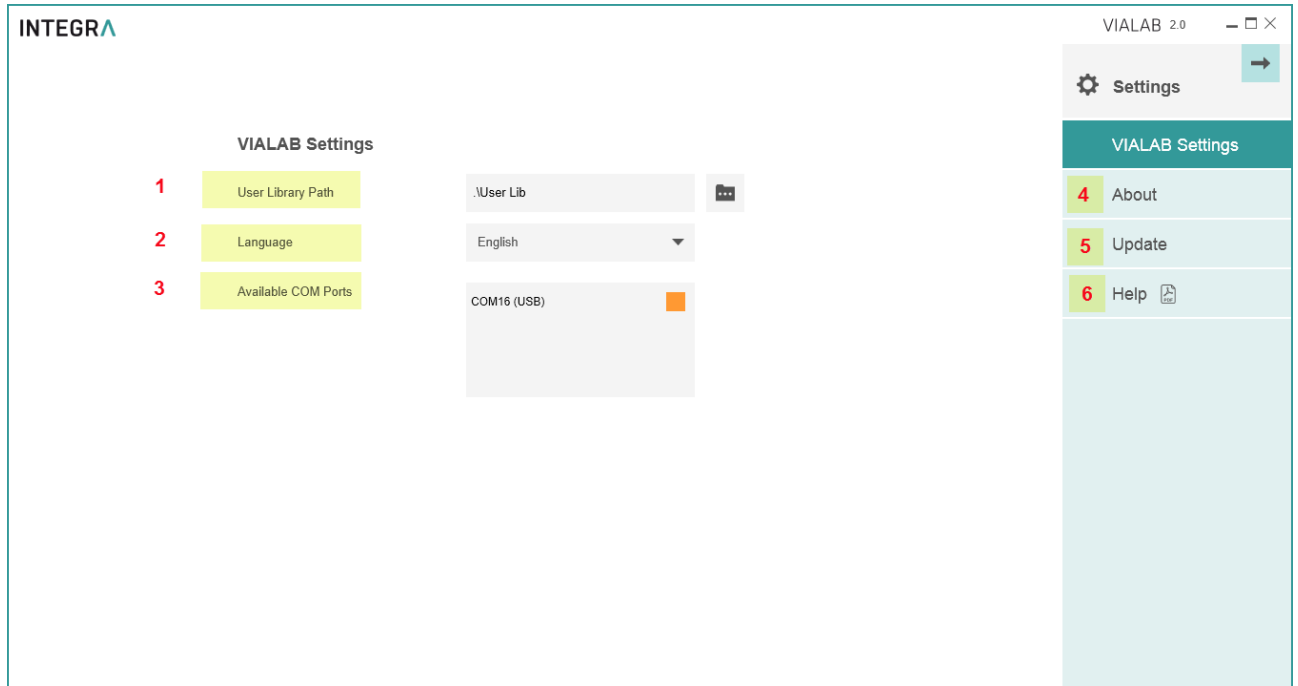

1 User Library Path: Define the storage path of your labware library.

2 Language: Choose your preferred language.

3 Available COM Ports: Shows all available communication ports.


4 About: Information about the VIALAB software version and the INTEGRA contact data

5 Update: Check for new updates, see 4.2.

6 Help: Download the latest VIALAB operating instructions.

4.2 Update VIALAB

If an update for VIALAB is available and you are connected to the internet you can check for new updates.

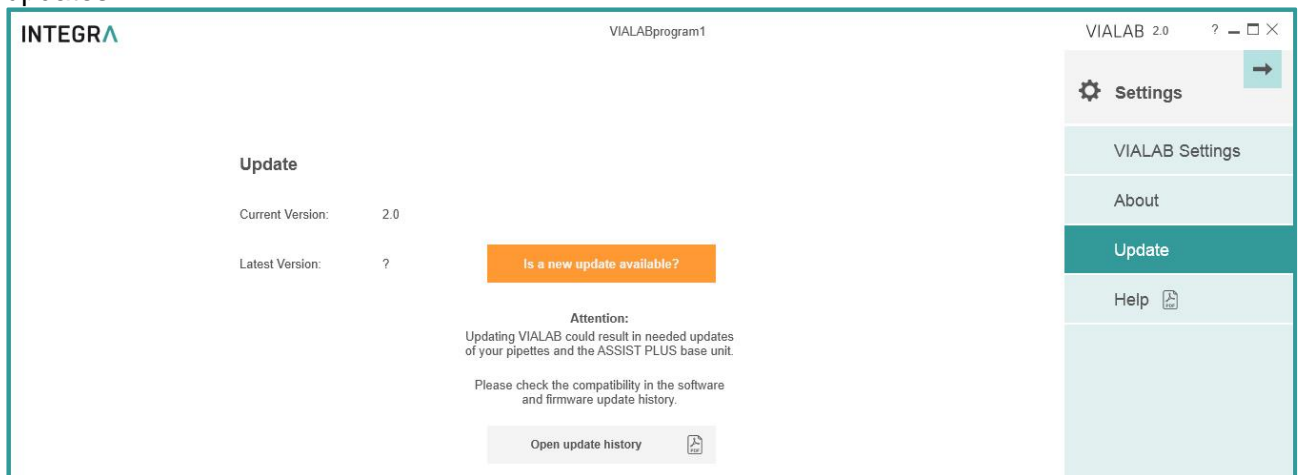

Click on Open update history to download information about software versions of VIALAB and VIALINK as well as firmware information of VIAFLO/VOYAGER pipettes and ASSIST / ASSIST PLUS.

5 Main Menu

5.1 Information and saving options



Click on the Menu button in the top left corner to open the main menu. Here information about the program is displayed and saving options are listed.

INTEGRA SerialDilution VIALAB 2.0

1 Information Insert additional information for this program.

2 Program Settings Program Name (on pipette) Description

3 Save SerialDilution 8 samples in 1.5 ml tubes are serially diluted with buffer in a 96 well plate.

4 Save As 14/15 Material: Ploette: VOYAGER - 8 ch - 300 ul 729/1000

5 Report Program overview:

Pipette	Tip	Deck Content
VOYAGER - 300ul - 8CH	300 ul GripTip, Sterile, Filter, Low retention	Initial Volumes Transfer

6 Close

7 New

8 Open

9 Sample Programs

10 Exit

- 1 Information:** change the Program name and add a description, if desired.
- 2 Program Settings:** define the clearance height and offsets of each deck position for your program.
- 3 Save:** save your program.
- 4 Save As:** save you program under a new name.
- 5 Report:** generate a printable PDF-report of your program, e. g. for standard operating procedures.
- 6 Close:** close your program.
- 7 New:** create a new program.
- 8 Open:** open an existing program, recently used programs are listed.
- 9 Sample Programs:** select and open a sample program.
- 10 Exit:** exit VIALAB

5.1.1 Program Settings


The following settings apply to each VIALAB program individually.

The screenshot shows the 'Settings for this program' window for 'VIALABprogram1'. The left sidebar contains navigation options: Information, Program Settings (selected), Save, Save As, Report, Close, New, Open, Sample Programs, and Exit. The main area displays settings for three different labware types: Multichannel Reservoir, Rack for 1.5 ml microcentrifuge tubes, and 96 well F-bottom. Each labware type has its own set of X, Y, and Z axis offset controls. At the bottom, there are three advanced menu options: Tip Ejection Type, Waste as target, and Copy Step - Height Adjustments. Red numbers 1 through 5 are placed next to the following settings: 1. Clearance Height (8.0 mm), 2. Deck position offsets (A, B, C), 3. Tip Ejection Type (Display/Hide), 4. Waste as target (Display/Hide), and 5. Copy Step - Height Adjustments (Automatic/Manual).

- 1 Clearance Height:** distance over the highest point of the labware at which the pipette moves over the labware
- 2 Deck Offsets:** deck positions A, B or C can be shifted in X, Y and Z axis for max. 2 mm
- 3 Tip Ejection Type:** choose whether to display or hide an additional option for the tip ejection. If activated in the tab "Tip Change" the tips can be ejected either at the ejector ledge (**Ledge**) or directly in the waste bin without touching the ejector ledge (**Waste**) during tip ejection.
- 4 Waste as target:** choose whether to display or hide the option to dispense directly into the waste bin. If activated the waste bin can be selected as target position in a Multiaspirate step or in Serial Dilution step as a waste position for the last aspiration.
- 5 Copy Step – Height Adjustments:** when copying steps the liquid level changes for the copied step. If tip travel is activated in a step, the heights are automatically adjusted relative to the new liquid level (Automatic). Alternatively, the heights of a copied step aren't changed and identical to the original step (Manual).

6 Programming


6.1 Overview

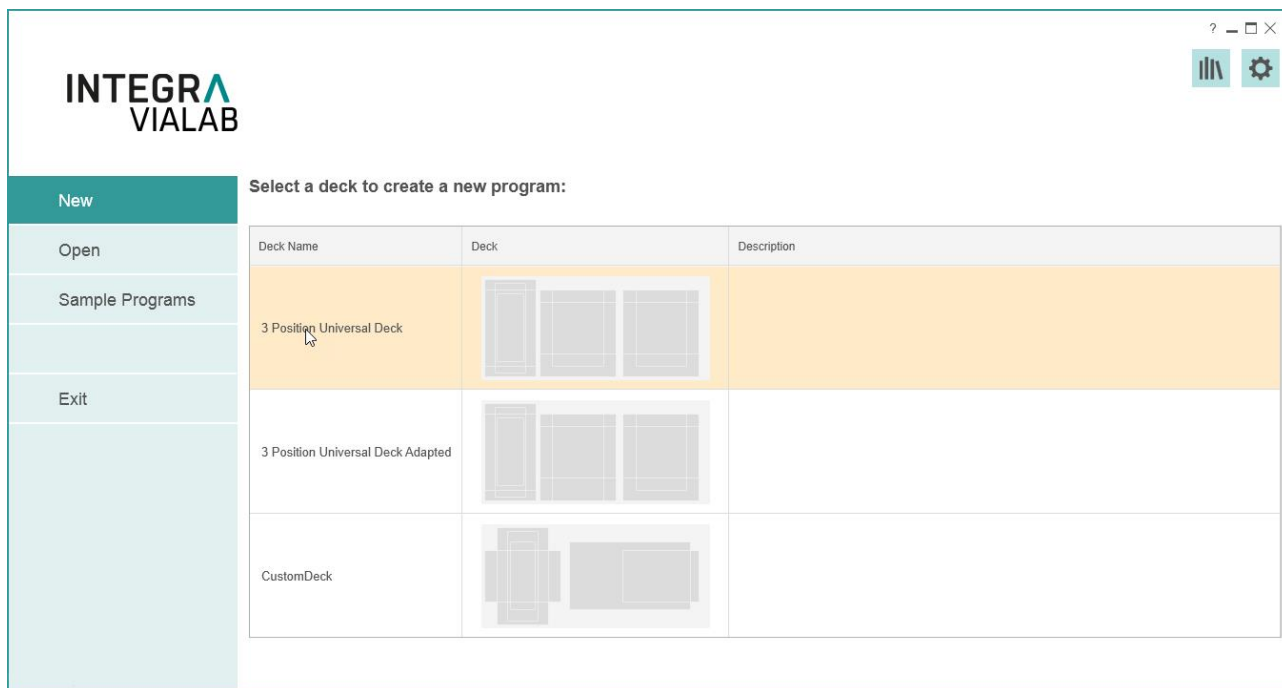

After starting VIALAB, you can create a new pipetting program, open an existing or recently used program or get started with a sample program.

To create a new program select first the deck you want to use, e.g. the 3 Position Universal Deck with one reservoir position and two spring loaded plate positions.

The program menu will be displayed:

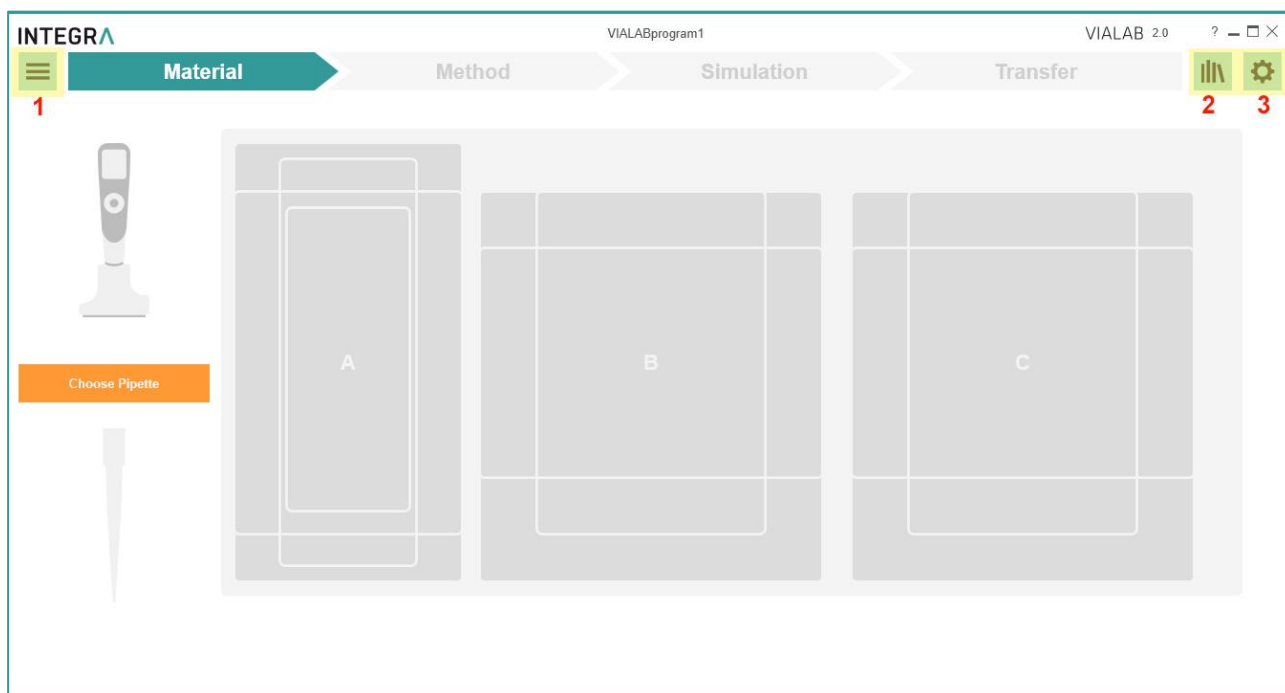


- 1 Main Menu:** Information and saving options of the defined program (see chapter 5)
- 2 Library:** Labware library listing available decks, pipettes, tips, reservoirs, plates, tube racks, and COMBI systems with COMBI bases which can be used for a program. You can also define new reservoirs, plates, tube racks and COMBI systems. Alternatively, import already defined labware files to the labware library.
- 3 Gear wheel:** Settings and updates of VIALAB.

The table below shows the different sections for creating a program. After one section is completed, the next can be defined.


Section	Description
Material	Choose the pipette with tips and the labware on the deck, e.g. the reservoir used as source container and the kind of well plates with orientation.
Method	Create your program that can contain up to 99 steps.
Simulation	Review a simulation of all steps in your pipetting program.
Transfer	Transfer the program to your pipette or reload an adapted program from the pipette.

6.2 Material

First click on the orange square to choose the pipette and afterwards select the GripTips. This enables the selection of the other labware used in the program.

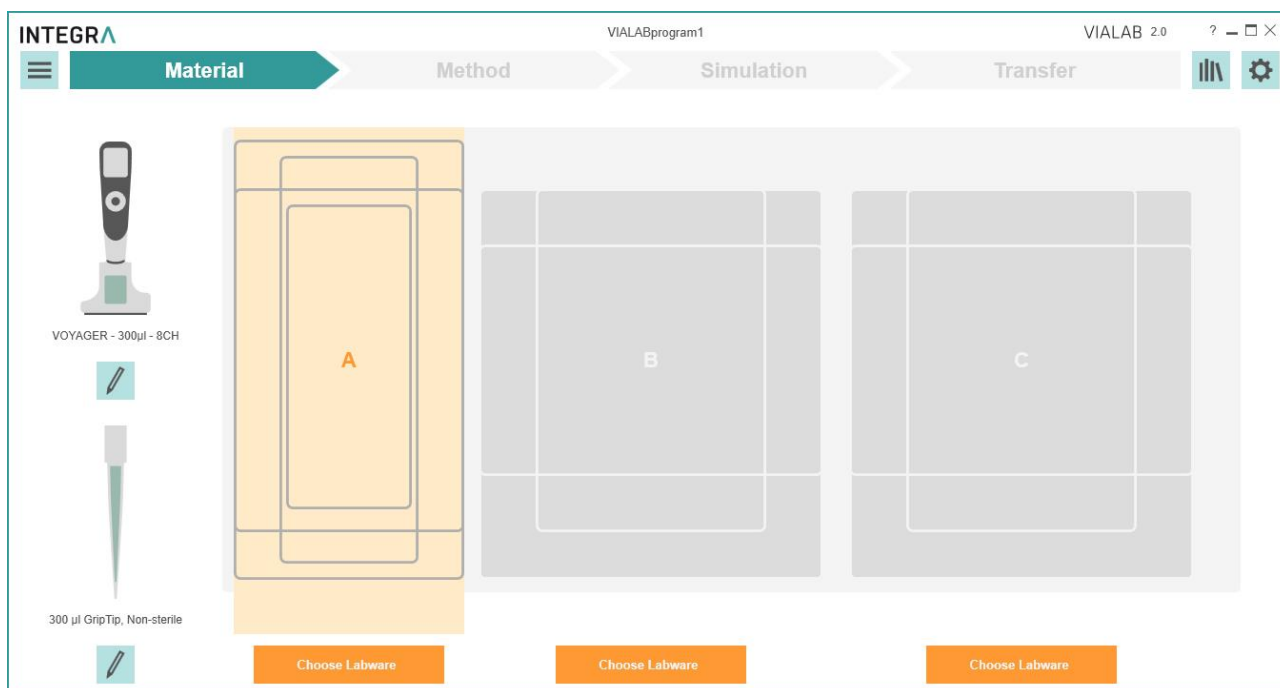

Click on an orange button to choose the predefined labware on all required positions.

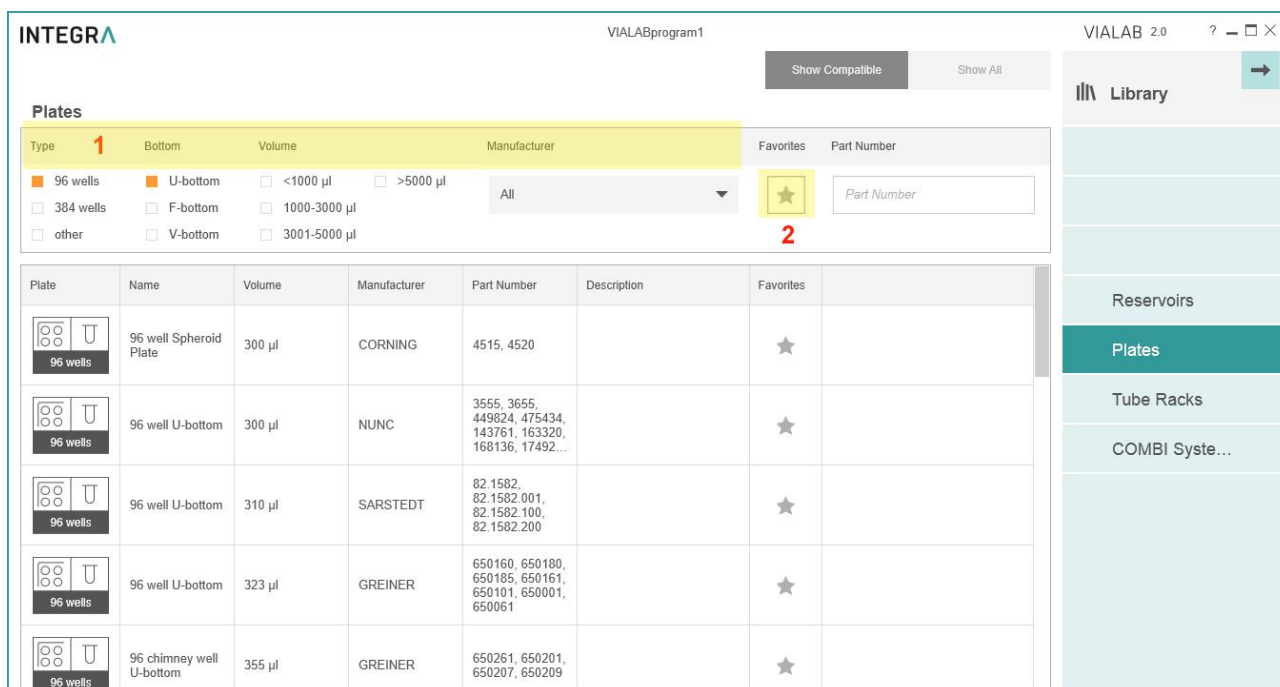


Select the desired properties or manufacturer (1) to quickly show reservoirs/plates/tube racks/COMBI Systems which are compatible with the selected pipette and deck position.

Click on the star (2) to show favorites only (to be defined in the labware library).

Scroll to one reservoir, plate, tube rack or COMBI System and select it by clicking on it.

Labware	Description
Reservoirs	Select an INTEGRA reservoir or define a new reservoir, see Library icon.
Plates	Select one of predefined, commercially available plates. Search for any plate or enter the part number directly, if known. Alternatively, define the dimensions of a new plate, see Library icon.
Tube Racks	Select a tube rack or define a new one, see Library icon.
COMBI System	To add a COMBI System, select a COMBI BASE and fill each selected column with modules, i.e. plates, tube racks or reservoirs and click on Apply COMBI System, see also section 7.5.

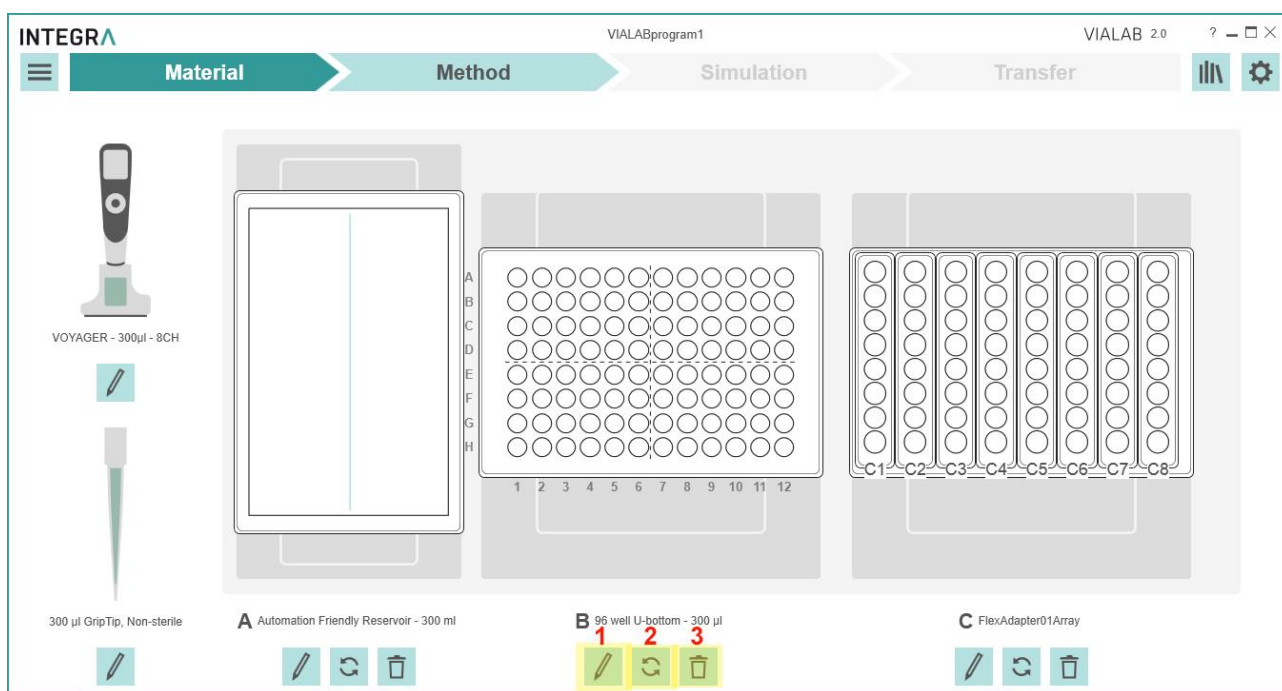


- 1 Edit:** Change the chosen labware.
- 2 Change orientation:** Rotate counterclockwise by 90° (all 4 possible orientations).
- 3 Delete:** Delete the selected labware on this position.

Note: Changing the material of an already defined program, e.g. pipette or plate orientation, will delete all subsequent program steps in the “Method” section, if they are not compatible with the previous settings.

6.3 Method

Your pipetting program is composed of several steps.

When you open the section Method the first time, the first step “Initial Volumes” is already preset and you are prompted to set the volumes of the labware which are already filled.

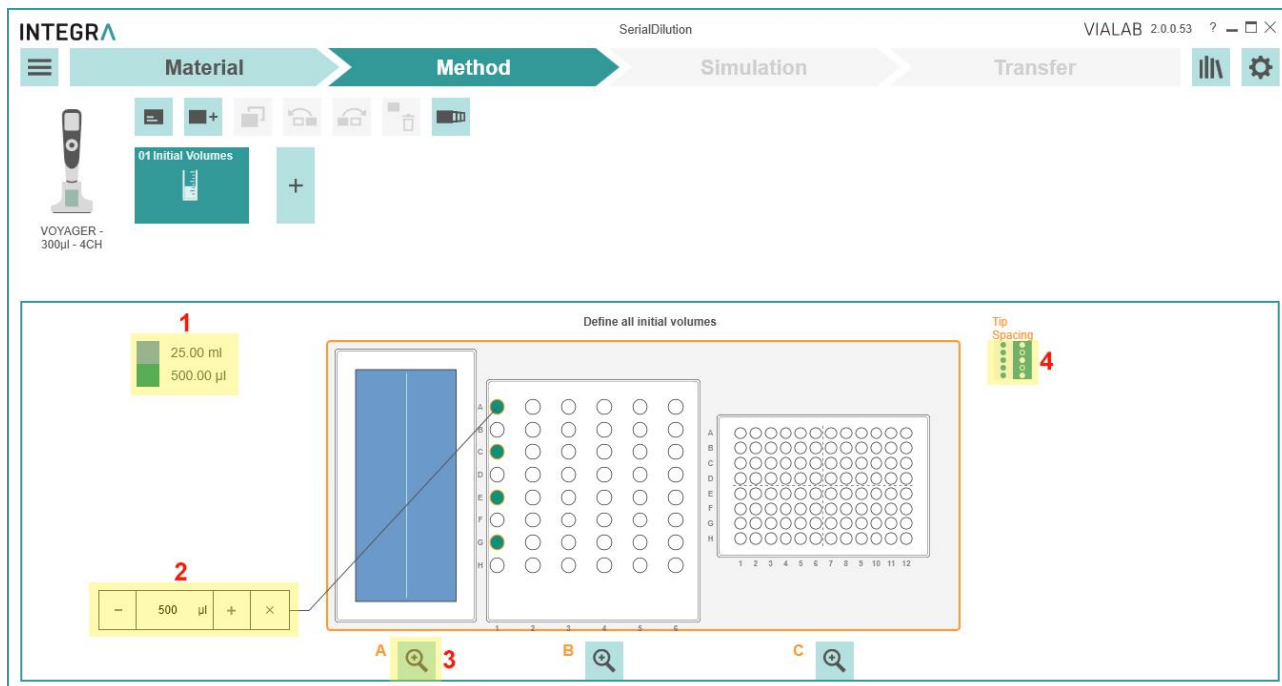


- 1 Color legend:** Shows the pre-fill volume of the labware on position A, B or C. Each color represents one pre-fill volume.
- 2 Pre-fill volume:** Click on the vessel which is already filled at the start of your program and enter the pre-fill volume.
- 3 Zoom:** Displays a close-up view of one deck position.
- 4 Tip spacing:** Choose a different tip spacing to change the filling pattern of your plate e.g. as shown with plate on position C (only for VOYAGER tip spacing pipettes).

Select the wells which are prefilled at the start of your program via click or drag and drop. In the pop-up window you type in the dedicated volume in the given unit to assign a volume.

To view or edit the assigned volume select wells and the current volume is displayed. If needed, edit the volume or delete it.



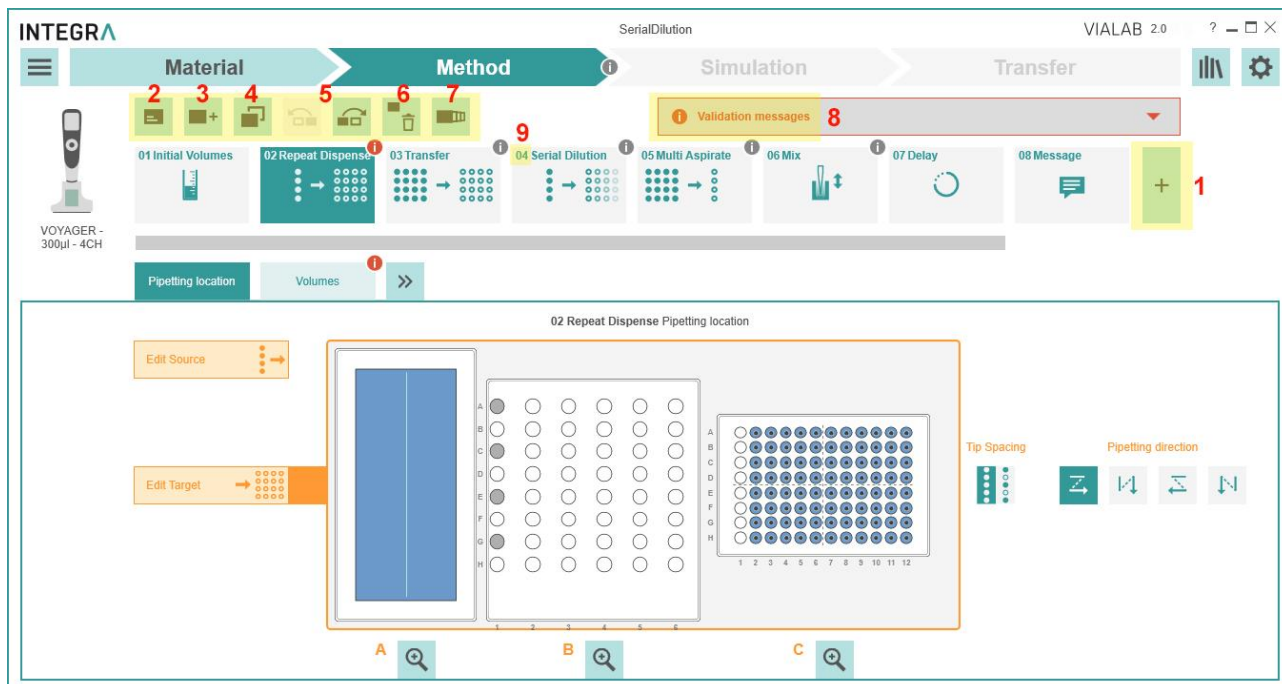

To add further steps click on a plus square (1, 3). To move, copy or delete steps click on the corresponding icon (4, 5, 6). The steps are processed sequentially according to the step number (9).




- 1 **Plus square:** add a new step to the end of the method.
- 2 **Add Description:** add a description for the current step. The first 15 characters are displayed at the bottom of the step icon.
- 3 **Add Step:** add a new step after a selected step to the method.
- 4 **Copy Step:** copy selected step to the end.
- 5 **Move left/right:** move the selected step forwards or backwards.
- 6 **Delete Step:** delete selected step.
- 7 **Multi-Selection:** click on this button to select multiple steps in parallel to copy or delete them.
- 8 **Validation messages:** click on this field for programming error details.
- 9 **Step number:** the step number indicates the step sequence in the program.


Following steps are available:

Symbol	Step	Description
	Transfer	Transfers an equal number of samples from a source to target (N:N), e.g. plate copy.
	Repeat Dispense	Dispenses from one source into multiple targets (1:N), e.g. reagent addition
	Serial Dilution	Aspirates sample from a source and performs a serial dilution in the defined number of wells.
	Multi Aspirate	Collects a number of samples from one source and transfers into one target (N:1).
	Mix	Mixes a sample by pipetting up and down in the selected vessel.


Symbol	Step	Description
	Delay	Defines a pause between the last and the next step in seconds.
	Message	Pauses the program and displays a message. Three lines with 12 characters each are available.
	Volume Change	Enter the filling volume of the vessels e.g. after manual filling steps within a method. The volume defined in this step is an absolute volume. It does not take into account previous liquid handling steps.
	AUX Port	Turning devices on and off which are connected to the auxiliary port of the ASSIST PLUS.



Within one step, adjust the following parameters:

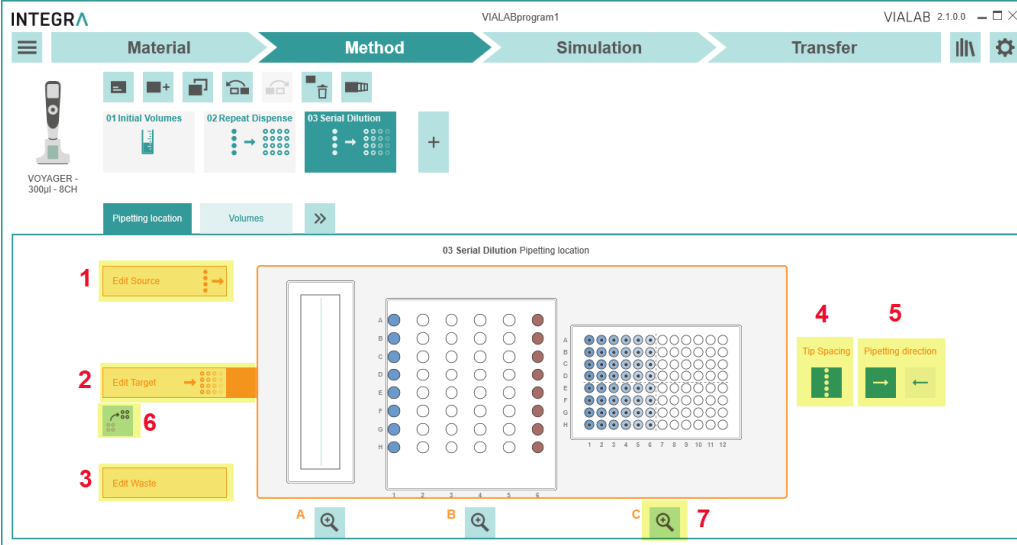Parameter	Description
Pipetting location	<ul style="list-style-type: none"> • Edit Source (1): Click on the button and then select the wells (for plates on upmost well of column) where the liquid should be aspirated from. It will be colored. • Edit Target (2): Click on the button and then select the wells (for plates on upmost well of column) where the liquid should be dispensed in. It will be marked with a black dot in the middle of the well. • Edit Waste (3): Only for step “Serial Dilution” and if under Volumes “Last aspiration back to “Waste” was selected. Click on the button to select waste vessel location.  <p>If needed, adapt the tip spacing (4) or the pipetting direction (5).</p> <p>For changing the selected wells deselect the chosen ones or move all selected wells together within one plate (6). For a detailed view of the plates, use the zoom-in function (7)</p>
Volumes	<p>Move the slider up or down (1), type in the volume (2) or use the “+” and “-” buttons to set the desired Volume(s) for the current step. An Air Gap can be set to aspirate a defined air volume after an aspiration of liquid.</p>

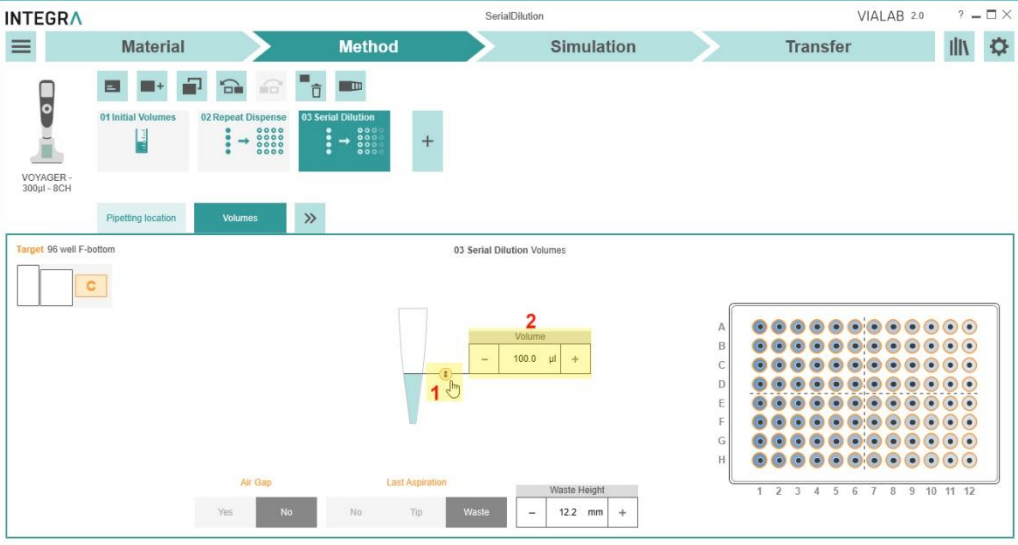Parameter	Description	
		
Volumes	Step Type	Volume description
	Transfer	Volume: Define whether the volume is variable or fix for all transfers. If variable, select each individual column and assign a volume.
	Repeat Dispense	<p>Post-Dispense Location: Select if the post-dispense should be purged in the source vessel of aspiration or should remain in the tips.</p> <p>Reuse Post-Dispense: Click on “Yes” if the post-dispense should remain in the tip and be used for the next aspiration.</p> <p>Dispense Type: Click on “Single” to aspirate and dispense the selected volume or on “Multi” to aspirate a multiple of the selected volume and perform several dispenses.</p> <p>Volume: Define whether the volume(s) to be dispensed (for every dispense step) is fix or variable. If variable, select each individual column and assign a volume. Set the pipetting direction. Additionally to the dispense volume set also the Pre-Dispense and the Post-Dispense volume.</p> <p>Note: We recommend a pre- and post-dispense volume of 8 % of the nominal volume of the pipette.</p>
	Serial Dilution	<p>Last Aspiration:</p> <ul style="list-style-type: none"> - “No”: the last aspiration is not performed, the liquid remains in the last row. - “Tip”: the last aspiration remains in the tips that are ejected into the Waste bin. - “Waste”: the last aspiration is purged in a dedicated liquid waste vessel. The location has to be defined in “Pipetting location”. <p>Volume: The sample volume that will be transferred from well to well.</p>
	Multi Aspirate	<p>Aspirate Type: Click on “Single” to aspirate only one volume or on “Multi” to perform several aspirations consecutively.</p> <p>Volume: Define whether the volume to be aspirated (for every aspirate step) is fix or variable. If variable, select each individual column and assign a volume. Set the aspiration volume.</p>
Speeds	<ul style="list-style-type: none"> • Aspiration and Dispense Speed: Set the aspiration and dispense speeds for this step. Set the slider to the appropriate speed from 1 (slow) to 10 (fast). 	

Parameter	Description
	<ul style="list-style-type: none"> • If necessary, adapt parameters for viscous and volatile liquids: Aspiration Delay: A delay in seconds after aspiration before moving out of the liquid. Dispense Delay: A delay in seconds after dispensing before moving out of the liquid. Exit Liquid Slowly: The tips move slowly out of the liquid (for aspiration and for dispense) • Movement Speeds: Changes the speeds of the Tower in X-direction, of the Pipette arm in Z-direction or the Pipette holder in Y-direction [1 (slow) to 10 (fast)]. • Pre-Wetting: If pre-wetting of the tips is desired activate it and set the number of cycles. All fresh tips are pre-wetting with the following aspiration volume and speeds.

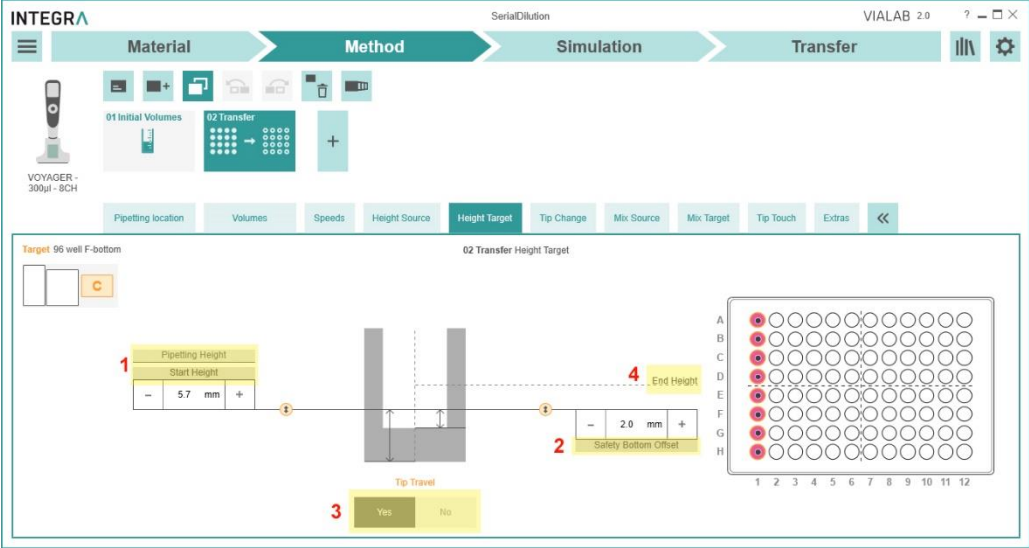Parameter	Description
Height Source Height Target	<p>The calculated liquid level in the well is displayed to ease height definition graphically.</p> <p>Click on Pipetting Height or Start Height (1) and move the line up or down to enter the heights for aspiration/dispense referring to the deck (deck height is 0 mm). The Safety Bottom Offset (2) defines the minimum distance of the tips to the bottom of the labware.</p> <p>If the Heights are variable, define the pipetting heights for every column.</p> <p>Tip Travel (3): Movement of tip during one aspiration/dispense to compensate the decrease/increase of the liquid level. Recommended for high pipetting volumes in vessels with small diameter.</p> <p>If Tip Travel is deactivated the aspiration/dispense is performed at a constant Pipetting Height (1).</p>  <p>If Tip Travel is activated aspiration/dispense starts at the Start Height (1) and ends at the End Height (4), which is calculated automatically according to the labware dimensions and pipetting volume.</p> <p>Note: It is possible to fine-tune the Pipetting Height / Start Height directly on the pipette to optimize the heights according to the used liquid (see ASSIST PLUS Operating Instructions chapter 5.4 “Adjust heights and positions”).</p>
Tip Change	<p>Set if and when tips should be changed.</p> <p>Continue next step without change: no tip ejection after program step. Same tips will be used for next step.</p> <p>After step complete: change tips after completing the current step</p> <p>After each transfer/dispense/dilution: change tips after each dispense</p> <p>Before each new aspiration: change tips before each new aspiration</p> <p>If the “Tip Ejection Type” is activated in the Program Settings (see chapter 5.1.1), the tips can be either positioned on the Ejector Ledge before ejection (Ledge) or are directly ejected into the middle of the waste bin (Waste). Tips may stack up and the waste bin needs to be emptied more often.</p>

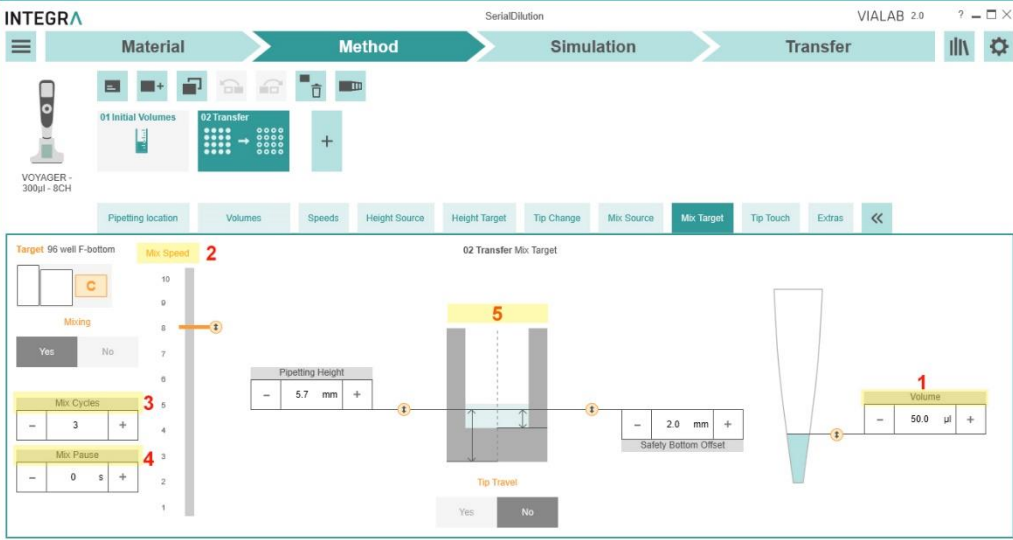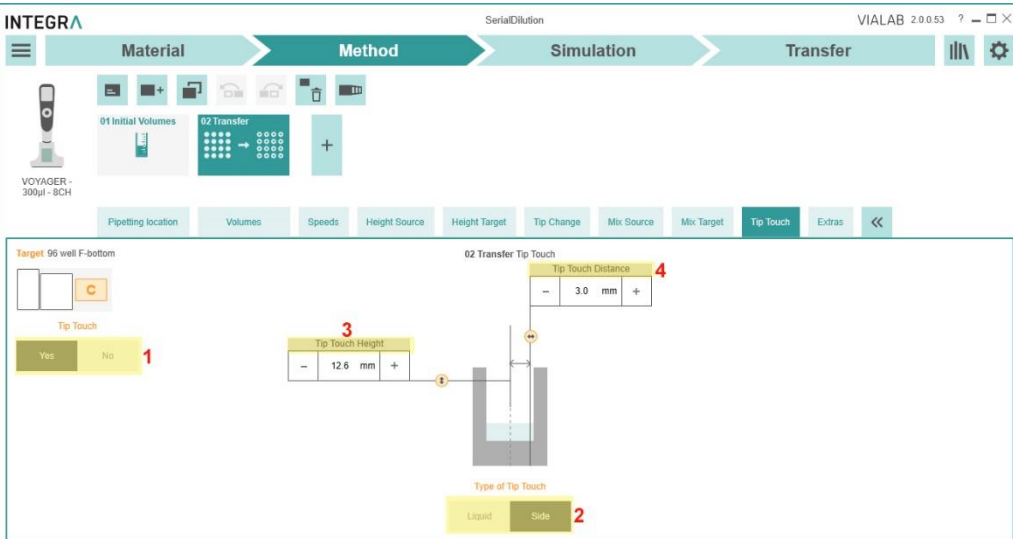Parameter	Description
Mix	<p>Define, if the sample should be mixed in the source well before aspiration and in the target well after dispense. Define the mix Volume (1), Speed (2), the number of Cycles (3) and if a Pause (4) should be performed in between mixing cycles. Define the mixing Heights (5).</p> 
Tip touch	<p>Set a tip touch (1) after dispensing to remove drops that may cling on the pipette tips. Set whether the tip touch should be performed into the Liquid or on the Side (2) of the wells. Enter the Tip Touch Height (3) and the Distance (4) to the vessel wall.</p> 

Parameter	Description
Extras	<div>Define a well Offset, if you want to pipette off the well center or even out of wells, e. g. for transwell plates or into the Slanted Plate Holder (PN 4510)</div> <div></div>


Note: If red exclamation marks appears, the program contains errors. Click on “Validation messages” for further details.

6.4 Simulation

In the simulation window the deck with the defined labware is displayed. On the top all steps are listed. Click on the play button to start the simulation of your program. Watch the whole simulation or navigate between steps or stop the simulation (1).

To speed up the simulation to 2-, 4- or 8 fold speed, select the appropriate button (2). An estimation of the Total Time and the Total Tip Consumption is listed (3).

The simulation is displayed as top (4) and side (5) view of the deck.

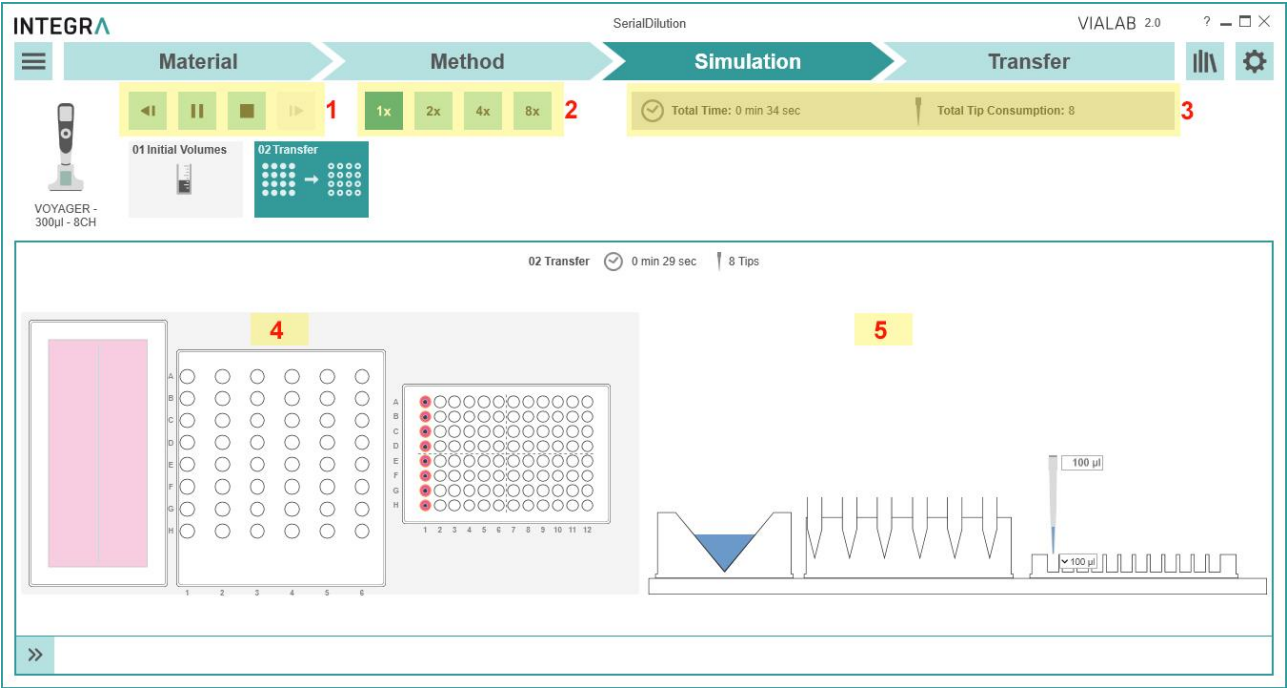

For all error-free programs a report is generated by clicking on the Main Menu and selecting "Report".
A report in pdf format with all settings of the program is generated.

6.5 Transfer

You should save your current program on your computer to have a backup or if you need to transfer the program to more than one pipette. Select the Main Menu button and click Save or Save As.

6.5.1 Connect pipette

Connection via Charging/Communication (USB) Stand (#4211)

Connect the USB cable to the stand and a free USB port on your computer. It is recommended to connect the stand to the mains power supply. Then hang the pipette onto the stand.

Go to the “Toolbox” of the pipette. In the menu “Communications” select “USB” on a VIAFLO/ VOYAGER electronic pipette. Stay in this menu for transfer.

Connection via Pipette Communication Module (#4221)


Pipette and PC need to be paired for the first connection. Please refer to MS Windows Help to learn how to pair Bluetooth devices. The following pairing code may be required: 12345.

Go to the “Toolbox” of the pipette. In the menu “Communications” select “ComModule”.

6.5.2 Transfer

In VIALAB switch to the section “Transfer” for saving the program to the pipette.

Upon establishing the connection successfully the pipette details, such as pipette model, serial number and firmware version, are displayed on the right side. Other VIALAB programs already saved on the pipette are listed.

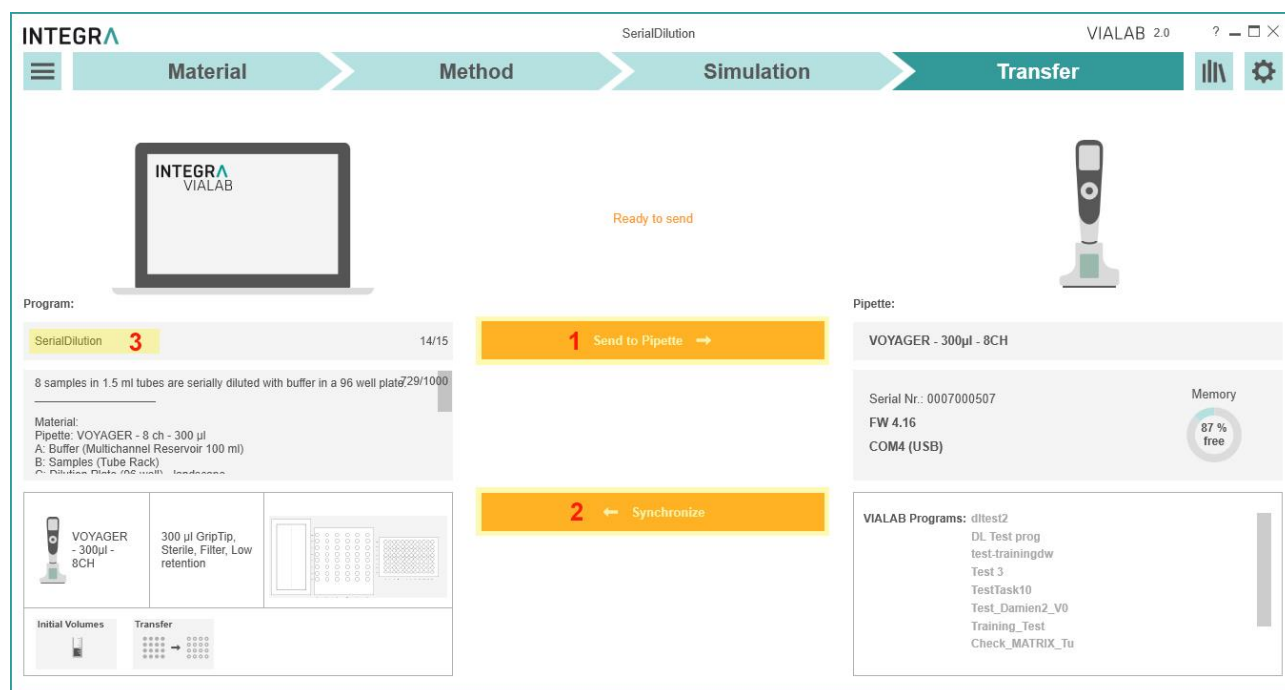

Send to Pipette (1) starts the transfer of the current program to the pipette. If certain parameters of the program are modified on the pipette (offsets, pipetting heights), synchronization of the program to VIALAB is possible (2).

To change the program display name on the pipette type in an alternative program name with max. 15 characters (3).

Successful saving or synchronization of the program to the pipette / to VIALAB is indicated by a temporary message at the bottom of the window.


7 Labware Library

You can access the labware library by clicking on the Library icon in the top right corner.

Within the labware library, you can define the dimensions of decks, reservoirs, plates, tube racks and COMBI Bases/Systems. Pipettes, GripTips and reservoirs from INTEGRA are already defined and cannot be modified.


7.1 Reservoirs

In the reservoir menu all predefined plates are listed. Select a desired property, e.g. Type, to quickly show already defined reservoirs of this property.

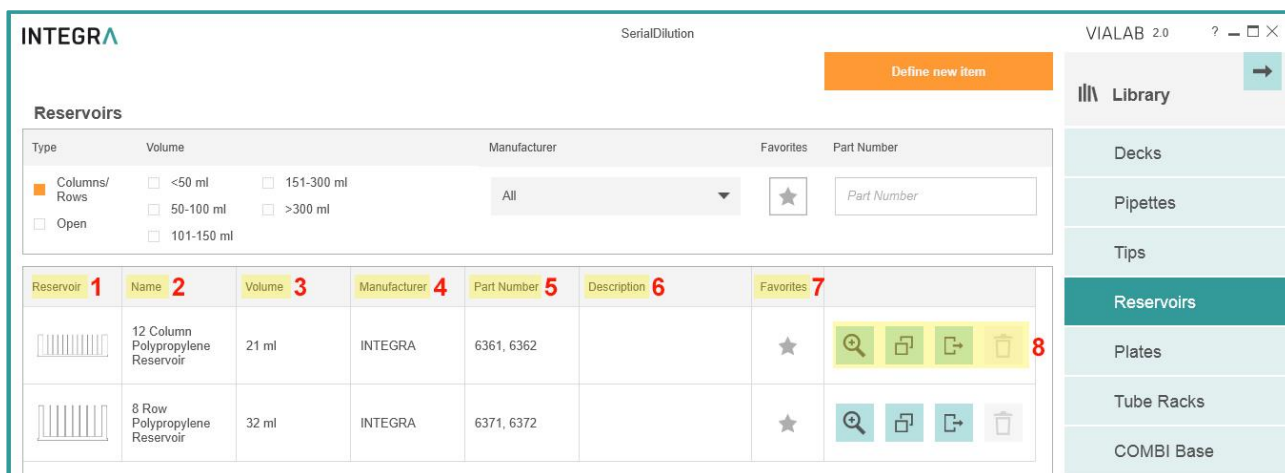


- 1 Reservoir:** the icon shows the number and shape of the wells and well bottom.
- 2 Name:** indicates the number and shape of the wells and well bottom, as well as other properties such as half-area.
- 3 Volume:** shows the maximal well volume of the reservoir or compartment as defined by the manufacturer.
- 4 Manufacturer:** name of the reservoir manufacturer
- 5 Part Number:** part number in accordance to the reservoir manufacturer
- 6 Description:** optional
- 7 Favorites:** select to mark as favorite. The material section in the program set-up has a filter to show only favorites.
- 8 Tools:** Edit/View details, Copy, Export or Delete the selected labware. It is not possible to edit or delete INTEGRA items. To modify other predefined labware, create a copy first and then edit the copy.

To add new reservoirs to the library, click on “Define new item”. In the Configure Reservoir menu enter a unique name, the manufacturer, the part number of the reservoir and a description, if applicable. The description can also be used for lab personal to identify the reservoir.

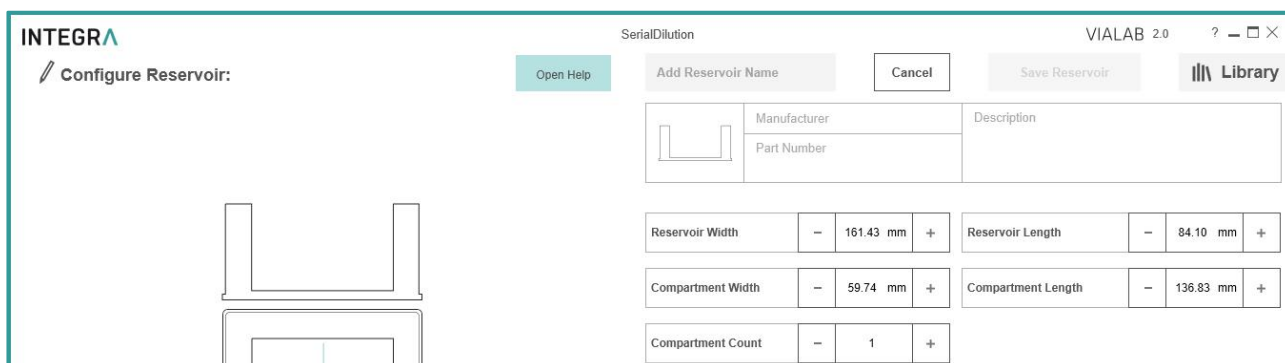

Enter the reservoir dimensions. You will be supported by live drawing of the defined reservoir.

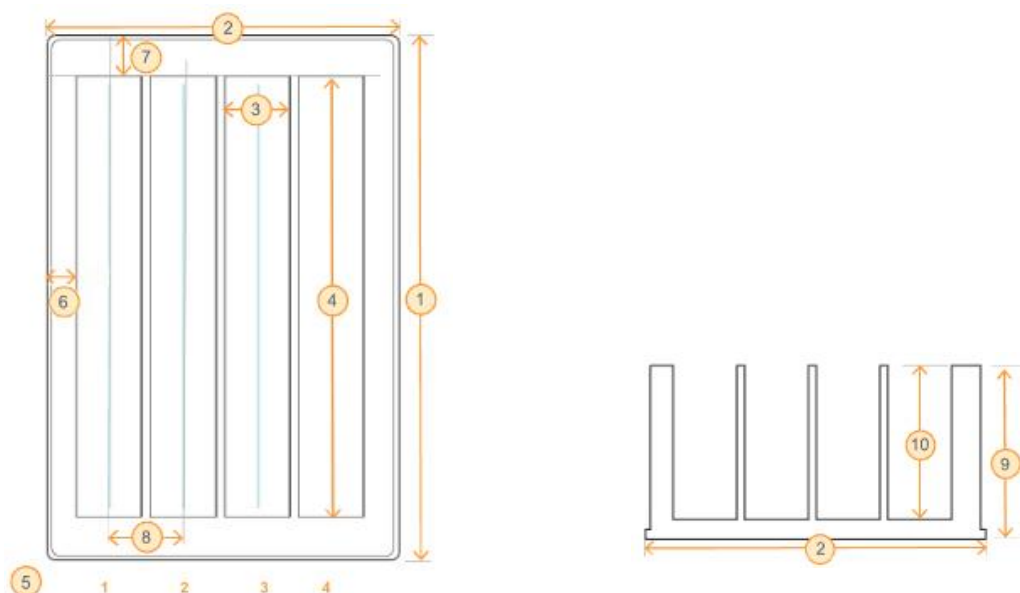

Note: Reservoirs are always set up with their compartment(s) in vertical direction.

Dimension	Description
Reservoir Width	The outside Width (1) of the reservoir including the bottom outside flanges, i.e. the displayed dimension from back to front.
Reservoir Length	The outside Length (2) of the reservoir including the bottom outside flanges, i.e. the displayed sideward dimension.
Compartment Width	The width of the compartments, i. e. the inner sideward dimensions (3).
Compartment Length	The length of the compartments, i. e. the inner distance from back to front (4).
Compartment Count	The number of the compartments (5), e. g. 4.
Left edge offset of First Compartment	The distance between the left outside edge of the reservoir and the left edge of the first compartment, including the bottom outside flange (6).
Top edge offset of First Compartment	The distance between the top outside edge of the reservoir and the top edge of the first compartment including the bottom outside flange (7).
Distance between Compartment Positions	The distance between the compartments, from the middle of one compartment to the middle of the next (8), e. g. 9 mm. This is the compartment width plus the spacing wall.
Type	Block. Rectangular flat bottom reservoirs can be defined.
Reservoir Height	The outside height (9) of the reservoir.
Compartment Depth	The inner depth (10) of the reservoir.
max. Volume	The maximal filling volume of one compartment.

7.2 Plates

In the plate menu all predefined plates are listed.

Type	Bottom	Volume	Manufacturer	Favorites	Part Number
<input type="checkbox"/> 96 wells	<input type="checkbox"/> U-bottom	<input type="checkbox"/> <1000 µl	<input type="checkbox"/> >5000 µl	All	★
<input type="checkbox"/> 384 wells	<input type="checkbox"/> F-bottom	<input type="checkbox"/> 1000-3000 µl			Part Number
<input type="checkbox"/> other	<input type="checkbox"/> V-bottom	<input type="checkbox"/> 3001-5000 µl			

Plate	Name	Volume	Manufacturer	Part Number	Description	Favorites
	12 well F-bottom	6325 µl	TPP	92012, 92412, 92112		★
	12 well F-bottom	6900 µl	CORNING	3336, 3512, 3513		★
	24 well F-bottom	1900 µl	CORNING	3337, 3524, 3526, 3527, 3473		★
	24 well F-bottom	3300 µl	GREINER	662102, 662160		★
	24 well F-bottom	3350 µl	TPP	92024, 92424, 92124		★

To add new plates to the library, click on “Define new item”.

In the Configure Plate menu enter a unique name (e.g. well number and shape), the manufacturer and the part number of the plate.

Configure Plate:

Open Help

Add Plate Name Cancel Save Plate Library

96 wells

Manufacturer Part Number Description

Row Count 8 Column Count 12

Well Shape - Top View Round Well Bottom V-bottom

Well Inner Diameter 8.00 mm

Center to Center Spacing - Rows 9.00 mm Center to Center Spacing - Columns 9.00 mm

Offset Row 1 11.24 mm Offset Column 1 14.38 mm

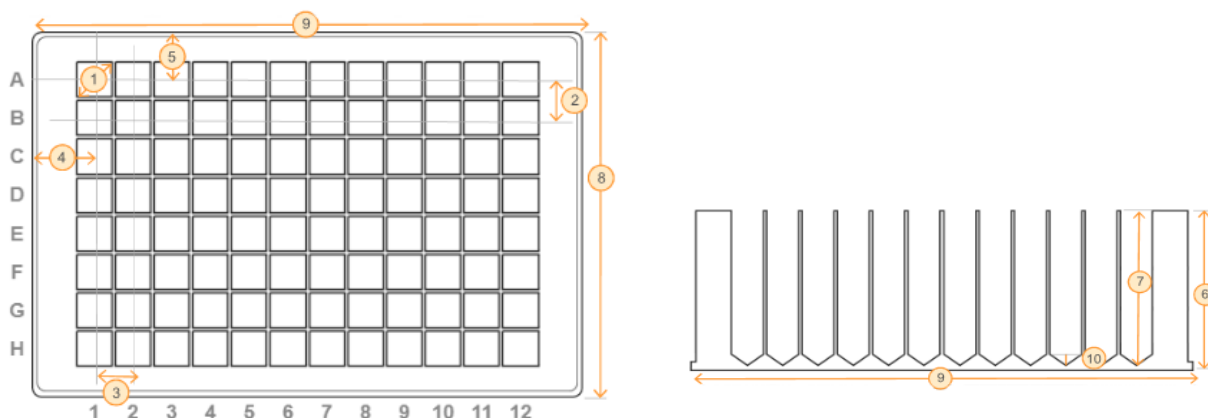
Plate Height 14.35 mm Well Depth 10.00 mm





Plate Width 85.48 mm Plate Length 127.76 mm

max. Well Volume 0 µl Well V-shape Depth 0.00 mm

Note: For plates in adapters, e. g. in PCR 96 or 384 well cooling plates, define the dimensions of the inserted plate in combination with the adapter as one plate.

Enter all plate dimensions. You will be supported by live drawing of the defined plate.



Dimension	Description	
Row Count	A 96 well microplate is arranged as 8 rows by 12 columns.	
Column Count	A 384 well microplate is arranged as 16 rows by 24 columns.	
Well Shape - Top View	Shape on the top of the well, is usually round.	
Well Bottom		U-bottom = round
		F-bottom = flat
		C-bottom = curved edges of flat bottom, e. g. for easy wash
		V-bottom = conical
Well Inner Diameter	Well diameter (1) on the top of the well	
Center to Center Spacing - Rows	The distance between the rows (2) /columns (3) from the middle of one well to the middle of the next well.	
Center to Center Spacing - Columns	96 well: ANSI/SLAS-standard ¹ : 9 mm 384 well: ANSI/SLAS-standard: 4.5 mm	
Offset Row 1	The distance between the top outside edge of the plate and the center of the first row of wells (4).	
	96 well	ANSI/SLAS-standard: 11.24 mm
	384 well	ANSI/SLAS-standard: 8.99 mm
Offset Column 1	The distance between the left outside edge of the plate and the center of the first column of wells (5).	
	96 well	ANSI/SLAS-standard: 14.38 mm
	384 well	ANSI/SLAS-standard: 12.13 mm
Plate Height	Plate height (6). ANSI/SLAS-standard: 14.35 mm ± 0.76 mm	
Well Depth	Depth in the center of the well (7).	
Plate Width	The shorter side (8), ANSI/SLAS-standard: 85.48 mm ± 0.5 mm	
Plate Length	The longer side (9), ANSI/SLAS-standard: 127.76 mm ± 0.5 mm	
max. Well Volume	The total, calculated well volume, not the working volume	
Well V-shape Depth	Depth from start of conical narrowing to well bottom only for V-bottom wells (10).	

¹ ANSI: American National Standards Institute, SLAS: Society for Laboratory Automation and Screening.

7.3 Tube Rack

To create new tube racks, select Tube Racks and click on “Define new item”. In the Configure Tube Rack menu enter a unique name, the manufacturer, the part number and a description, if applicable.

INTEGRA

Configure Tube Rack:

Open Help

Add Tube Rack Name

Cancel

Save Tube Rack

Library


96 wells





Manufacturer

Part Number

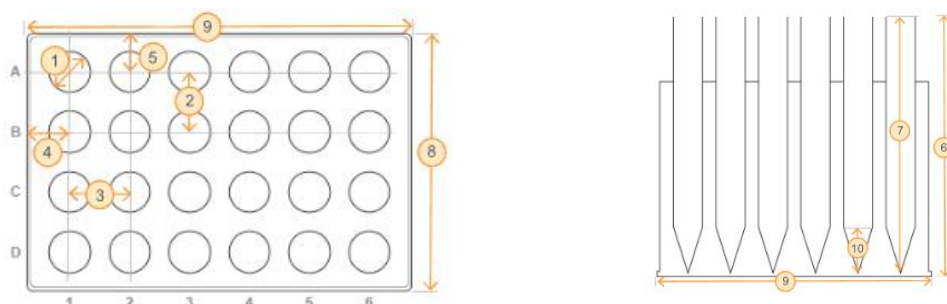Description

Enter all tube rack dimensions and the tube type. You will be supported by live drawing of the defined rack.




Dimension	Description	
Row Count	Number of rows.	
Column Count	Number of columns.	
Tube Shape - Topview	Shape on the top of the tubes, is usually round.	
Tube Shape - Bottom		U-bottom = round
		F-bottom = flat
		C-bottom = curved edges of flat bottom
		V-bottom = conical, including typical conical tube form with round bottom
Tube Inner Diameter	Tube inner diameter in the middle of the tube (1).	
Center to Center Spacing - Rows	The distance between the rows, from the middle of one tube to the middle of the next tube (2).	
Center to Center Spacing - Columns	The distance between the columns, from the middle of one tube to the middle of the next tube (3).	
Offset Column 1	The distance between the left outside edge of the rack and the center of the first column of tubes (4).	
Offset Row 1	The distance between the top outside edge of the rack and the center of the first row of tubes (5).	
Tube Rack Height	Total height of tube rack and tubes (6).	
Tube Depth	Depth in the center of the tube (7).	
Tube Rack Width	Width of the rack(8). ANSI/SLAS-standard: 85.48 mm ± 0.5 mm.	
Tube Rack Length	Length of the rack (9). ANSI/SLAS-standard: 127.76 mm ± 0.5 mm.	
Max. Tube Volume	The total, calculated tube volume, not the working volume.	
Tube V-shape Depth	Depth from start of conical narrowing to tube bottom, e. g. 3.1, only for V-bottom wells (10).	

7.4 COMBI Base

A COMBI Base is a platform where a module, i.e. a reservoir, plate or tube rack, can be placed on. It only acts as an adapter and you cannot choose a COMBI Base alone as a pipetting location (only a COMBI System see chapter 7.5)

To create a new COMBI Base, select COMBI Base and click on “Define new item”. In the Configure COMBI Base menu enter a unique name, the manufacturer, the part number of the COMBI Base and a description, if applicable.

Enter all COMBI base dimensions. You will be supported by live drawing of the defined COMBI base.

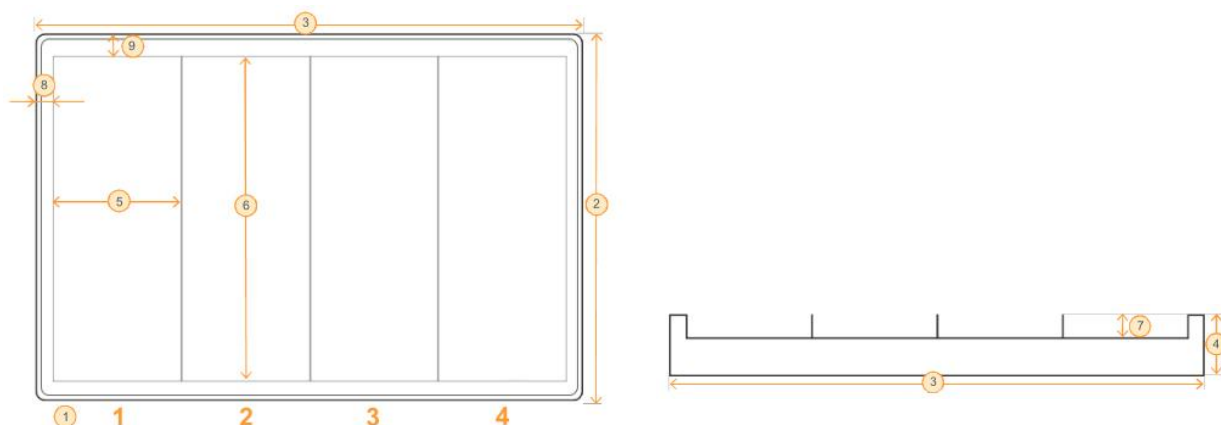

Dimension	Description
Division	Regular: a repeated number of identical sections are lined up. Irregular: sections are positioned freely within the outer dimensions of the base.
Sections	Number of sections (1).
Base Width	Width of the base (2). ANSI/SLAS-standard: 85.48 mm ± 0.5 mm.
Base Length	Length of the base (3). ANSI/SLAS-standard: 127.76 mm ± 0.5 mm.
Base Height	Total height of the base (4).
Sections Width	The width of a section (5).
Section Length	The length of a section (6).
Section Height	The height of a section (7). Negative values indicate that the bottom of the section is deeper than the whole COMBI base.
Left edge offset of section	The distance between the left outside edge of the base and the left edge of the first section, including the bottom outside flange (8).
Top edge offset of section	The distance between the top outside edge of the base and the top edge of the first section including the bottom outside flange (9).

Offset Row 1 The distance between the top outside edge of the rack and the center of the first row of tubes (5).

7.5 COMBI System

A COMBI System consists of a COMBI Base and individual modules, i.e. a reservoir, plate or tube rack, which are placed on the COMBI Base.

To create a new COMBI System, click on “Define new item”. In the Combine COMBI System menu select a COMBI Base.

Enter a unique name, the manufacturer, the part number and a description of the COMBI System, if applicable.

COMBI System: Select a module

Open Help

Slider Adapter_System

Cancel

Save COMBI System

Library

Manufacturer

Part Number

Adapter for up to 6 sliders (Base from PN 4540 and PN 4541)

☐ Plates

☐ Tube Racks

☐ Reservoirs

All

Part Number

	Name	Manufacturer	Part Number	Description	Favorites...
	Slider 0.5 ml microcentrifuge tubes	INTEGRA		Slider for 0.5 ml microcentrifuge tubes (see PN 4541)	★
	Slider 1.5 ml microcentrifuge tubes	INTEGRA		Slider for 1.5 ml microcentrifuge tubes (see PN 4540)	★
	Slider 2 ml microcentrifuge tubes	INTEGRA		Slider for 2 ml microcentrifuge tubes (see PN 4540)	★

1 Name: Unique name under which the COMBI system is stored.

2 Save COMBI system: Save button.

Click in a section and select a Module. Continue with the next section until all sections are filled. You will be supported by live drawing of the defined COMBI System.

Click on Save COMBI System.

7.5.1 Use of a Slanted Plate Holder (PN 4510) in VIALAB

The Slanted Plate Holder can accommodate plates in three different angles: 10°, 20° and 30°.


To assemble a Slanted Plate Holder with a plate select in “Material” deck position **B** or **C** on the 3 Position Universal deck. Choose the menu option “COMBI” on the right menu bar and select the COMBI Base “Slanted Plate Holder” at the desired angle.

Select either the landscape or portrait deck section on the Slanted Plate Holder and choose a compatible plate from the list on the right side.

Choose “Apply COMBI System” and use it for your program.

If you want to reuse a COMBI System, you can also create a new COMBI System in the labware library and save it with an appropriate name. Then the already assembled COMBI System (adapter and plate) can be used in your method.

Note: A change of the angle needs a change of the COMBI system used in the Material section. All steps in the Method section will be deleted then. We advise to define the appropriate angle in a test program first, before programming a lot of steps, which might be lost.

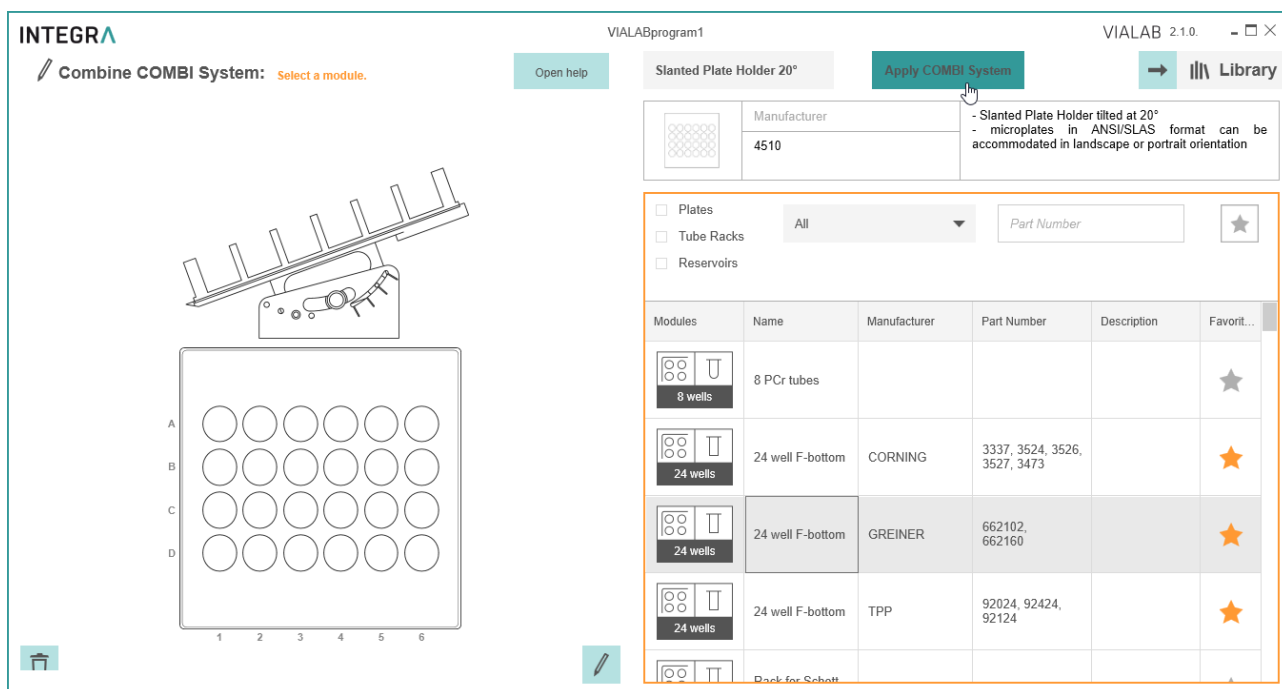

In the labware library, the Slanted Plate Holder (20°) is assembled e.g. with a 24 well plate.

Height Adjustments for the Slanted Plate Holder

When using the Slanted Plate Holder the pipetting (and mixing) heights are different between the individual columns.

The default heights can be adjusted by selecting each column individually and adjusting the height.

The lowest selectable height in the displayed graphic in refers to the lowest accessible height in the tilted well in one column.

To access the inner rim and the very bottom of a well it might be needed to set an X-Offset ("Extras" tab) to shift the position of the tip within the well.

Depending on the tilt angle and orientation of the adapter, the plate type, the tip size and position you want to access in the plate the offset is different.

8 Trouble Shooting

Topic	Issue	Remedy
Labware Library	I cannot place my labware item on the deck, although it is defined correctly in the labware library.	<ul style="list-style-type: none"> • Make sure the selected labware item is compatible with your pipette.
		<ul style="list-style-type: none"> • Make sure your labware item fits to the selected deck position e.g. a tube rack cannot be placed on deck position A.
Labware Library / Method	I created my labware item and placed it on the deck, but I cannot fill it with volume.	<ul style="list-style-type: none"> • Make sure the Max. Well Volume is > 0 µl.
Connection between pipette and computer	I cannot get a connection between pipette and VIALAB.	<ul style="list-style-type: none"> • Make sure the pipette is in communication mode (see chapter 6.5.1.) • Make sure the Programming Stand or device is connected to the computer. • ComModule: The Pipette Communication Module and PC need to be paired if connected for the first time. Refer to Windows Help to learn how to pair Bluetooth devices. The passcode 12345 may be required. • Make sure the drivers for the Programming Stand are installed. While the Programming Stand is connected, go to your System Settings -> Device Manager. Check under Ports (COM & LPT) if there is an entry "USB Serial Port (COM X)". If there is not, install the drivers manually.
Synchronization of programs from pipette to VIALAB.	The "Synchronize" button is not activated.	<ul style="list-style-type: none"> • Make sure that no changes were made in VIALAB after saving the program to the pipette.
Pipette status in section Transfer.	I connected the pipette and its status is "Pipette is incompatible".	<ul style="list-style-type: none"> • Make sure that the correct pipette model matching the current program is connected. You can only transfer programs to the pipette model used in the program.
Program Transfer To Pipette	The pipette is connected, after "Send To Pipette" the program is not saved on the pipette.	<ul style="list-style-type: none"> • Make sure, that you do not use any special characters in your program name on the pipette (Transfer Screen).