

MEDIAJET Petri Dish Compatibility List

DOCUMENT NAME: OI_MJ_Petri_dish_compatibility_list_V10.docx

Document Change History

Vers.	Status:	Name:	Date:	Change:
V08	Released	ULe	01.11.2018	Discontinued plates removed, new tables, file renamed.
V09	Released	ULe	05.04.2019	Dishes from FL Medical and Kord Valmark added.
V10	Released	ULe	20.12.2019	Dishes from LabServ, Corning added. New format.

Table of Contents

1	Introduction	1
2	General Recommendations	2
3	Petri dishes Ø90/Ø100 - MEDIAJET standard/vario	3
4	Petri dishes Ø55/Ø60 - MEDIAJET vario	5
5	Petri dishes Ø35 - MEDIAJET vario	6

1 Introduction

This document is intended to support sales representatives and customers working with MEDIAJET devices and applications.

For historical reasons, almost all PDs are not specifically designed for automated use. They are neither standardized. You can find an extensive variation of different Ø90/Ø100, Ø55/Ø60 or Ø35 dishes on the market. The actually measured dimension of the Petri dishes deviate from the nominal dimensions given by the manufacturer. We found variations of heights and diameters within the same product ID, even depending on the production location, and within the same lot.

Thus a general compatibility of a locally available Petri dish with MEDIAJET or MEDIAJET vario cannot be guaranteed by INTEGRA Biosciences.

For an exact specification of the dimensions for Petri dishes compatible with MEDIAJET please refer to the latest Operating Instructions (OI), 6.2 Petri dish specifications. In addition, this document lists a set of Petri dishes already successfully tested for mechanical processability and used on MEDIAJET by INTEGRA Biosciences.


As manufactures can change the dimensions of their Petri dishes without prior informing INTEGRA Biosciences, this list can only serve as an indication and is subjected to continuous changes.



We recommend to test the application with the desired dishes and original media in the customer's lab.

INTEGRA

2 General Recommendations

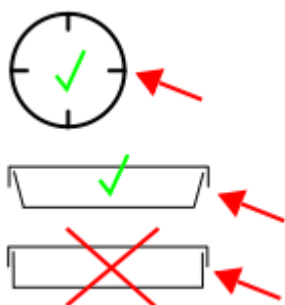


Use dishes with ventilation camps on the inside of the dish lid. Ventilation camps prevent the formation vacuum, the base can fall down.


Dishes with stacking ring can be used, because the gripper mechanism gently holds the plates without shearing them away.

The base of the plate must be conical.


The following dimensions are critical for an optimal performance of the MEDIAJET:

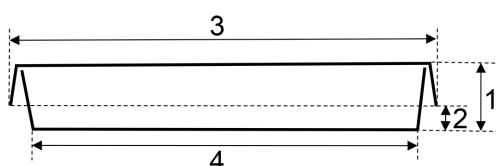

1. **Dish height:** 12-25 mm
2. **Ground clearance:**
 ≥ 6 mm (for \varnothing 90 mm dishes) or
 ≥ 5 mm (for \varnothing 60 and \varnothing 35 mm dishes)
3. **Lid outer diameter**
4. **Base outer diameter**

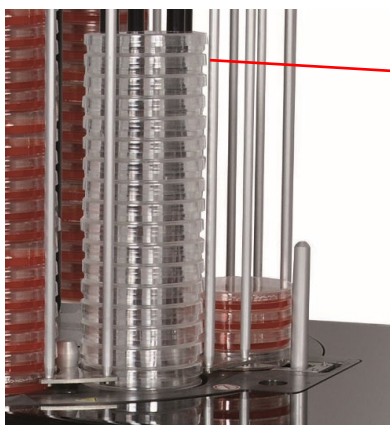

If the lid outer diameter (at bottom) is too large, the lid get stuck between the carousel rods.

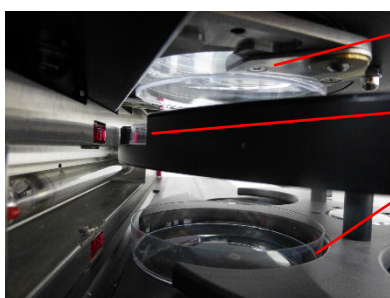

If the ground clearance is too low, the gripper cannot grip the base without lid.

If the lid outer diameter is too small, the lid falls through the upper rotor position and the dish is not opened.

If the base outer diameter is too large, the base does not lie flat in the rotor and dispensing is not possible without spilling.

Petri dish manufacturer nominal size	Lid outer diameter (3) (mm)	Base outer diameter (4) (mm)	Carousel	Rotor	Conversion set
90 – 100 mm	93.2 – 94.6	≤ 88.4	90L	90L	90L
	91.7 – 93.1		90	90 (standard or biplate)	90 (required for MJ vario only)
	90.7 – 92.1	≤ 87.6		90S	
55 – 60 mm	58.0 – 59.4	≤ 54.8	60	60	60
	56.6 – 58.0	≤ 51.2		60S	60S
35 mm	37.8 – 39.2	≤ 35.4	35	35	35

Special configurations for other Petri dish formats are available on request.

3 Petri dishes Ø90/Ø100 - MEDIAJET standard/vario

3.1 Petri dishes for 90 configuration (MEDIAJET standard PN 103005/103006 or MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
BD Falcon (Corning)	351029	100 x 15 mm		for Ø90 mm: Part No. 103 020 or 103 021	for Ø90 mm: Part No. 103 271 or biplate Part No. 113 806	for Ø90 mm, MEDIAJET vario only: Part. No. 113 051
Bibby-Sterilin	101VR20	90 mm				
Corning Gosselin	SH93-101 BH93B-102 (sterile)	90 x 14.2 mm				
	SB93-101 SB93-102 (sterile) BP93-101 (bag) BP93-102 (sterile, bag)	100 x 15 mm	90S rotor is also suitable			
	BT93-02	100 x 25 mm	deep			
	Fisherbrand	FB0875712 standard	100 x 15 mm			
FL Medical, Italy	29052, 29062 (sterile)	90 x 16.2 mm				
Greiner bio-one	633180 633102 (heavy) 633161 (heavy, sterile)	94 x 16 mm				
	664102 (heavy) 664160 (cellstar, ster.)	100 x 20 mm				
Kord Valmark, USA	2909	100 x 25 mm	deep			
Nunc	4031	100 x 25 mm	deep			
Phoenix Biomedical	018	90 x 15 mm				
Sarstedt AG	82.1473	92 x 16 mm				
	82.1472	92 x 16 mm				
VWR	391-0468, 391-0608 391-0469 (sterile)	90 x 16.2 mm				

Note: Product numbers in brackets are not tested, but the dimensions are identical (e.g. sterile and non-sterile version).

3.2 Petri dishes for 90S configuration (MEDIAJET standard PN 103005/103006 or MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
BIO-BIK	I-90	90 mm		for Ø90 mm: Part No. 103 020 or 103 021	for Ø90S mm: Part No. 103 272	for Ø90 mm, MEDIAJET vario only: Part. No. 113 051
LabServ, New Zealand	LBS60018X	90 x 20 mm				
Sarstedt AG	82.1473	92 x 16 mm				
	82.1472	92 x 16 mm				

3.3 Petri dishes for 90L configuration (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
				for Ø90L mm Part No. 113 022	for Ø90L mm Part No. 113 460	for Ø90L mm Part. No. 113 055

3.4 Petri dishes for biplate (two-compartment) rotor configuration (MEDIAJET vario PN 113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
Bibby-Sterilin	502V	90 mm		for Ø90 mm Part No. 103 020 or 103 021	for Ø90 mm biplate Part No. 113 806	for Ø90 mm Part. No. 113 051
Corning Gosselin	BI93D-101	100 x 15 mm				
FL Medical, Italy	29057 (29067)	90 x 16.2 mm				
Greiner bio-one	635102 635161 (sterile)	94 x 16 mm				
Sarstedt AG	82.1195	92 x 16 mm				
VWR	391-0492	90 x 16.2 mm				

4 Petri dishes Ø55/Ø60 - MEDIAJET vario

4.1 Petri dishes for 60 configuration (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
BD Falcon (Corning)	351007	60 x 15 mm		for Ø60 mm Part No. 113 020	for Ø60 mm Part No. 113 271	for Ø60 mm Part. No. 113 052 Conversion set
Corning-Gosselin	BP53-06 BP53-03 (sterile)	55 x 14.2 mm	Base of 54.1 mm too large for 60S rotor. BP50: same dimensions, no vents.			
Greiner bio-one	628102 628160 (cellstar, ster.) 628161 (sterile)	60 x 15 mm				
Kord Valmark, USA	2916	60 x 15 mm				
Sarstedt AG	82.1194	60 x 15 mm	Set stack height to 16.1 instead of 15.8			
VWR	391-3655	60 x 15 mm				

4.2 Petri dishes for 60S configuration (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
				for Ø60 mm Part No. 113 020	for Ø60S mm Part No. 113 276	for Ø60 mm Part. No. 113 056

5 Petri dishes Ø35 - MEDIAJET vario

5.1 (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
Greiner bio-one	627102 627160 (cellstar, ster.) 627161 (sterile)	35 x 15 mm		for Ø35 mm	for Ø35 mm	for Ø35 mm Part. No. 113 053
Nunc	153 066	35 x 10 mm		Part No. 113 021	Part No. 113 272	