

## MEDIAJET Petri Dish Compatibility List

DOCUMENT NAME: OI\_MJ\_Petri\_dish\_compatibility\_list\_V21.docx

### Document Change History

Vers.	Status:	Name:	Date:	Change:
V19	Released	GPo	09.01.2023	Tarsons 90 mm and 60 mm dishes added
V20	Released	GPo	25.09.2023	VWR 100 x 15 mm #25384-342 added BIOFIL MCD100090 added GREINER #664161 (VWR PN 82050-918) added VWR PN 25384-302 added
V21	Draft	GPo	01.05.2024	Kord Valmark #2901 60x15 added Fisher FB0875713 and FB08757150 Bi-Plate added NUNC #150288 (60x15 mm) added Fisher FB0875711 added

### Table of Contents

1	Introduction	1
2	General Recommendations	2
3	Petri dishes Ø90/Ø100 - MEDIAJET standard/vario	4
4	Petri dishes Ø55/Ø60 - MEDIAJET vario	10
5	Petri dishes Ø35 - MEDIAJET vario	12

## 1 Introduction

This document is intended to support sales representatives and customers working with MEDIAJET devices and applications.

For historical reasons, almost all PDs are not specifically designed for automated use. They are neither standardized. You can find an extensive variation of different Ø90/Ø100, Ø55/Ø60 or Ø35 dishes on the market. The actually measured dimension of the Petri dishes deviate from the nominal dimensions given by the manufacturer. We found variations of heights and diameters within the same product ID, even depending on the production location, and within the same lot.

Thus a general compatibility of a locally available Petri dish with MEDIAJET or MEDIAJET *vario* cannot be guaranteed by INTEGRA Biosciences.


For an exact specification of the dimensions for Petri dishes compatible with MEDIAJET please refer to the latest Operating Instructions (OI), 6.2 Petri dish specifications. In addition, this document lists a set of Petri dishes already successfully tested for mechanical processability and used on MEDIAJET by INTEGRA Biosciences.

As manufactures can change the dimensions of their Petri dishes without prior informing INTEGRA Biosciences, this list can only serve as an indication and is subjected to continuous changes.



We recommend to test the application with the desired dishes and original media in the customer's lab.

## 2 General Recommendations

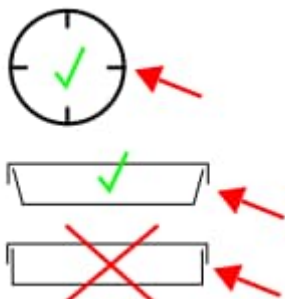


Use dishes with ventilation camps on the inside of the dish lid. Ventilation camps prevent the formation vacuum, the base can fall down.


Dishes with stacking ring can be used, because the gripper mechanism gently holds the plates without shearing them away.

The base of the plate must be conical.


The following dimensions are critical for an optimal performance of the MEDIAJET:

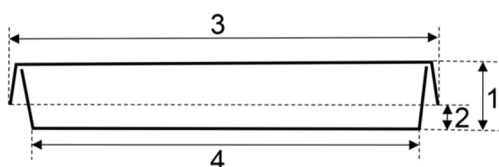

1. **Dish height:** 12-25 mm
2. **Ground clearance:**  
 ≥ 6 mm (for Ø 90 mm dishes) or  
 ≥ 5 mm (for Ø 60 and Ø 35 mm dishes)
3. **Lid outer diameter**
4. **Base outer diameter**

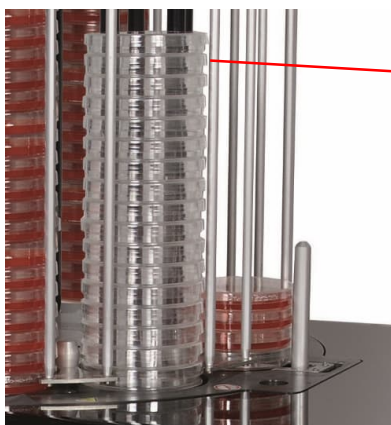

If the lid outer diameter (at bottom) is too large, the lid gets stuck between the carousel rods.

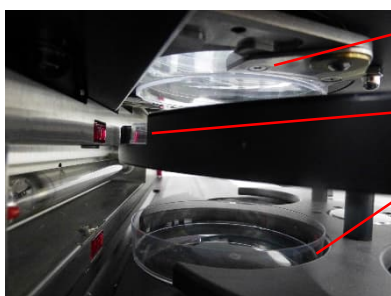

If the ground clearance is too low, the gripper cannot grip the base without lid.

If the lid outer diameter is too small, the lid falls through the upper rotor position and the dish is not opened.

If the base outer diameter is too large, the base does not lie flat in the rotor and dispensing is not possible without spilling.

Petri dish manufacturer nominal size	Lid outer diameter (3) (mm)	Base outer diameter (4) (mm)	Carousel	Rotor	Conversion set
90 – 100 mm	93.2 – 94.6	≤ 88.4	90L	90L	90L
	91.7 – 93.1		90	90 (standard or biplate)	
	90.7 – 92.1	≤ 87.6	90S	90S	90 (required for MJ vario only)
55 – 60 mm	58.0 – 59.4	≤ 54.8	60	60	60

# INTEGRA

	56.6 – 58.0	≤ 51.2		60S	60S
35 mm	37.8 – 39.2	≤ 35.4	35	35	35

Special configurations for other Petri dish formats are available on request.

### 3 Petri dishes Ø90/Ø100 - MEDIAJET standard/vario

#### 3.1 Petri dishes for 90 configuration (MEDIAJET standard PN 103005/103006 or MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
BD Falcon (Corning)	351029	100 x 15 mm		for Ø90 mm: Part No. 103 020 or 103 021	for Ø90 mm: Part No. 103 271  or <b>biplate</b> Part No. 113 806	for Ø90 mm, MEDIAJET vario only: Part. No. 113 051
Bibby-Sterilin	101VR20	90 mm				
Corning Gosselin	SH93-101 BH93B-102 (sterile)	90 x 14.2 mm				
	SB93-101 SB93-102 (sterile) BP93-101 (bag) BP93-102 (sterile, bag)	100 x 15 mm	90S rotor is also suitable			
	BT93-02	100 x 25 mm	deep			
	Fisherbrand	FB0875712 standard	100 x 15 mm			
FL Medical, Italy	29052, 29062 (sterile)	90 x 16.2 mm				
Greiner bio-one	633180 633102 (heavy) 633161 (heavy, sterile) 632181	94 x 16 mm	Tested by IUS			
	664102 (heavy) 664160 (cellstar, ster.)	100 x 20 mm				
Kord Valmark, USA	2909	100 x 25 mm	deep			
Nunc	4031	100 x 25 mm	deep			
Phoenix Biomedical	018	90 x 15 mm				
Sarstedt AG	82.1473	92 x 16 mm				
	82.1472	92 x 16 mm				
VWR	391-0468, 391-0608 391-0469 (sterile)	90 x 16.2 mm				
VWR	89107-632	100 x 25 mm	Deep dishes, tested by IUS, Spacer			

# INTEGRA

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
			under filling nozzle bracket required.			
VWR	25384-342	100 x 15 mm	Tested by IUS, plate runs smoothly			
VWR	25384-302	100 x 15 mm	Teach in successful only with UV enabled. 155mm 10 dish height.			
Techno Plas (AU)	10603001	90 x 14 mm	Ground clearance (6 mm) just within specifications.	for Ø90 mm: Part No. 103 020 or 103 021	for Ø90 mm: Part No. 103 271  or <b>biplate</b> Part No. 113 806	for Ø90 mm, MEDIAJET vario only: Part. No. 113 051
Techno Plas (AU)	TECTPG9014	90 x 14 mm				
NEST	753001	100 x15 mm	Tested by IUS, works well, 10 dish height: 157 mm			
Tarsons	460095	90 x 14 mm	Runs smoothly			

**Note:** Product numbers in brackets are not tested, but the dimensions are identical (e.g. sterile and non-sterile version).

## 3.2 Petri dishes for 90S configuration (MEDIAJET standard PN 103005/103006 or MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
BIO-BIK	I-90	90 mm		for Ø90S mm: Part No. 103 022 or 103 023	for Ø90S mm: Part No. 103 272	for Ø90 mm, MEDIAJET vario only: Part. No. 113 051
LabServ, New Zealand	LBS60018X	90 x 20 mm				
Sarstedt AG	82.1473	92 x 16 mm				
	82.1472	92 x 16 mm				
EIKEN Chemical (JP)	AU3110	87.2 x 14.25 mm	Lid diameter (90.7 mm) just within 90S rotor specifications. Centering unit (PN 103048) for planar stacking of PDs is strongly recommended			
EIKEN Chemical (JP)	AU3240	87.3 x 14 mm	Lid diameter (90.4 mm) is below 90S rotor specifications. PDs are still properly processed. Centering unit (PN 103048) for planar stacking of PDs is strongly recommended			
Sansei (JP)	01-002 CL	90 x 15 mm	Lid diameter (90.4 mm) is below 90S rotor specifications. PDs are still properly processed. Centering unit (PN 103048) for planar stacking of PDs is strongly recommended			
Sansei (JP)	01-014 SCL	90 x 15 mm	Lid diameter (90.5 mm) is below 90S rotor specifications. PDs are still properly processed. Centering unit (PN 103048) for planar stacking of PDs is strongly recommended			
Sansei (JP)	01-004 FX	90 x 20 mm	Deep PDs. Lid diameter (90.5 mm) is below 90S rotor specifications. PDs are still properly processed. Centering unit (PN 103048) for planar stacking of PDs is strongly recommended			
BIOFIL	MCD100090	90 x 15 mm	PDs run smoothly with 90S configuration. Centering unit 103048 to optimize planar stacking is recommended			

# INTEGRA

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
GREINER	664161 (VWR PN 82050-918)	100 x 20 mm	210mm 10 dish height Dishes run error-free	for Ø90S mm: Part No. 103 022 or 103 023	for Ø90S mm: Part No. 103 272	for Ø90 mm, MEDIAJET vario only: Part. No. 113 051
Fisher	FB0875713	100 x 15 mm	Tested at IUS Runs well - Requires 90S components			
Fisher	FB0875711	100 x 25 mm	Tested at IUS Runs well - Requires 90S components 10 dish height 243 mm			

### 3.3 Petri dishes for 90L configuration (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
				for Ø90L mm Part No. 113 022	for Ø90L mm Part No. 113 460	for Ø90L mm Part. No. 113 055



## 3.4 Biplate (two-compartment) Petri dishes (for MEDIAJET vario with biplate option PN 113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
Bibby-Sterilin	502V	90 mm		for Ø90 mm Part No. 103 020 or 103 021	for Ø90 mm <b>biplate</b> Part No. 113 806	for Ø90 mm  Part. No. 113 051
Corning Gosselin	BI93D-101	100 x 15 mm				
FL Medical, Italy	29057 (29067)	90 x 16.2 mm				
Greiner bio-one	635102 635161 (sterile)	94 x 16 mm				
Sarstedt AG	82.1195	92 x 16 mm				
VWR	391-0492	90 x 16.2 mm		for Ø90S mm: Part No. 103 022 or 103 023	for Ø90S mm: Part No. 103 272	for Ø90 mm  Part. No. 113 051
Sansei (Japan)	01-006 H	90 x 15	Lid diameter (90.5 mm) is below 90S rotor specifications. PDs are still properly processed. Centering unit (PN 103048) for planar stacking of PDs is strongly recommended			
Fisher	FB08757150 Bi-Plate	100x15	Tested at IUS Runs well – requires 90S components (10 dish height = 148 mm)			

## 4 Petri dishes Ø55/Ø60 - MEDIAJET vario

### 4.1 Petri dishes for 60 configuration (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
BD Falcon (Corning)	351007	60 x 15 mm		for Ø60 mm Part No. 113 020	for Ø60 mm Part No. 113 271	for Ø60 mm Part. No. 113 052
Corning-Gosselin	BP53-06 BP53-03 (sterile)	55 x 14.2 mm	Base of 54.1 mm too large for 60S rotor. BP50: same dimensions, no vents.			
Greiner bio-one	628102 628160 (cellstar, ster.) 628161 (sterile)	60 x 15 mm				
Kord Valmark, USA	2916	60 x 15 mm				
Sarstedt AG	82.1194	60 x 15 mm	Set stack height to 16.1 instead of 15.8			
VWR	391-3655	60 x 15 mm				
FisherBrand	FB0875713A	60 x 15 mm	Tested by IUS. Plate runs well.			
Weber Scientific	3051-15	60 x 15 mm	Tested by IUS. Plate runs well.			
Tarsons	460061	60 x 15 mm	Tested at ICH. Plate runs smoothly			
Kord Valmark	2901	60 x15 mm	Tested at IUS 10 dish height: 147 mm			
NUNC	150288	60 x 15 mm	Plates run smoothly.			

### 4.2 Petri dishes for 60S configuration (MEDIAJET vario PN 113000/113001/113002)

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
Tritech Research	T3308	60 x 15 mm	Tested by IUS, plate runs well, no issues	for Ø60 mm	for Ø60S mm	for Ø60 mm

# INTEGRA

				Part No. 113 020	Part No. 113 276	Part. No. 113 056
--	--	--	--	---------------------	---------------------	----------------------

# INTEGRA

## 5 Petri dishes Ø35 - MEDIAJET vario

### 5.1 MEDIAJET vario PN 113000/113001/113002

Manufacturer	Product no./type	Nominal Size	Remark	Carousel	Rotor	Conversion set
Greiner bio-one	627102 627160 (cellstar, ster.) 627161 (sterile)	35 x 15 mm		for Ø35 mm	for Ø35 mm	for Ø35 mm Part. No. 113 053
Thermo Scientific (Nunc)	153 066	35 x 10 mm		Part No. 113 021	Part No. 113 272	