

INTEGRA



PIPETBOY comfort Accurate pipetting at your fingertips

PIPETBOY comfort

Flexible, lightweight and always ready to use

Charger and stand – all in one

Either with or without a pipette, on the table or on the wall – the utmost flexibility and ease of use for your laboratory. Your PIPETBOY comfort is on stand-by and fully charged, whenever you need it.



Light as a feather

A lightweight 190 grams, the strain is taken out of even prolonged pipetting work. Hands, arms and shoulders remain fatigue-free, enhancing your comfort.



Ergonomically designed

The PIPETBOY comfort sits snugly in the hand naturally matching its shape. The rounded handle also fits perfectly in small hands. Pipetting has never been so comfortable.



Flexible and compatible

The silicone Pipette Mount guarantees a sure grip on practically all available plastic or glass pipettes in the 1 – 100 ml volume range.



Maximum Safety

The easy-to-change sterile filter provides your PIPETBOY comfort with reliable protection from liquid contamination and damage. And if you do happen to aspirate critical liquids with the pipette, the PIPETBOY comfort silicone Pipette Mount, the Pipettor Nose Piece and the filter rubber are all autoclavable while the sterile filter can be replaced. So easy and so safe.




Battery indicator

An illuminated yellow LED shows when the battery needs to be charged. While the battery is being charged, the LED flashes to show that the contact has been correctly established. You know the status of your rechargeable battery at all times. When the PIPETBOY comfort is connected directly to the mains adapter, you can still pipette while it is charging. The PIPETBOY comfort is always on standby!



From rapid pipetting to sensitive precision dosing

You can control the speed of pipetting continuously via the knurled screw. Precise dosing is also possible by applying different amounts of pressure to the triggers.

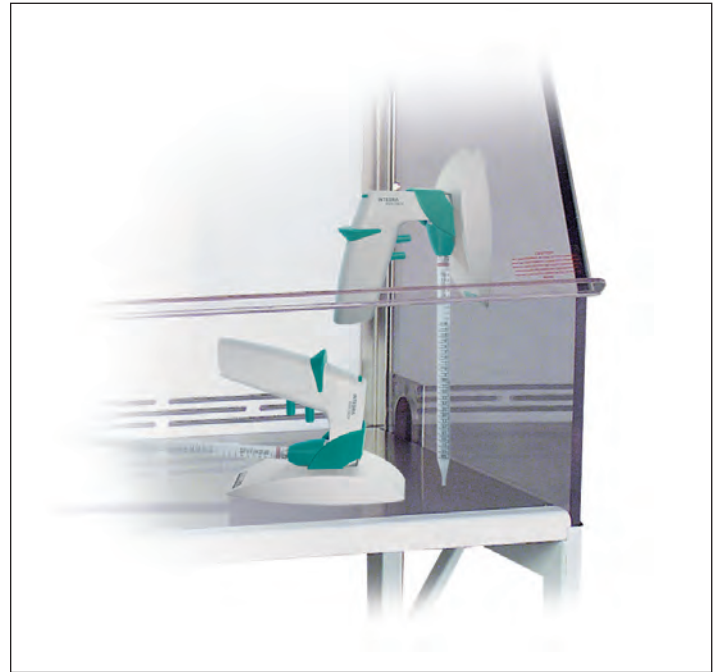

Held perfectly in all positions

Whether on the wall or tabletop, inside or outside a hood, the PIPETBOY comfort is held perfectly – with or without a pipette.

Wall mounting:

Magnet, adhesive strip or screw – the charger can be easily and securely mounted on any wall. The aerodynamic shape of the charger ensures minimal air turbulence in the laminar flow and reduces the risk of contamination.

Tabletop mounting:

Simply place the charger stand on top of your table – 4 rubber feet ensure stability.



Powerful

2 NiMH rechargeable batteries deliver optimum performance ensuring maximum power and durability with minimal weight. Fully charged batteries provide 8 hours of continuous operation. The PIPETBOY comfort can either be connected directly to the mains adapter or recharged via the charging console.

PIPETBOY comfort

Ordering information


Description	Part No.
PIPETBOY comfort classic , incl. charging station, adapter, sterile filter, charger mounting kit and two batteries	156 000
PIPETBOY comfort volcano , as Art. No. 156 000	156 001
PIPETBOY comfort pacific , as Art. No. 156 000	156 002
PIPETBOY comfort sahara , as Art. No. 156 000	156 003
PIPETBOY comfort royal , as Art. No. 156 000	156 004

Ordering information for single and spare parts

No.	Description	Part No.
1	Filter rubber	153 225
2	Filter sterile, 0.45 µm	153 015
	Filter unsterile, 0.45 µm	153 016
	Filter sterile, 0.2 µm	153 608
	Filter unsterile, 0.2 µm	153 607
3	Pipette Mount, silicone	151 020
4	Pipettor Nose Piece (classic)	156 100
	Pipettor Nose Piece (volcano)	156 101
	Pipettor Nose Piece (pacific)	156 102
	Pipettor Nose Piece (sahara)	156 103
	Pipettor Nose Piece (royal)	156 104
5	Battery-Set, 2 Pieces (1.2 V AAA NiMH/1000 mAh)	156 151
6	Battery cover Set (classic)	156 110
	Battery cover Set (volcano)	156 111
	Battery cover Set (pacific)	156 112
	Battery cover Set (sahara)	156 113
	Battery cover Set (royal)	156 114
7	Charging station (classic)	156 050
	Charging station (volcano)	156 051
	Charging station (pacific)	156 052
	Charging station (sahara)	156 053
	Charging station (royal)	156 054
8	Charger mounting kit (tape, magnet, base pads, screw)	156 150
9	Adapter:	
	US version	156 630
	(100 – 240 VAC, 50/60 Hz)	
	EU version	156 631
	UK version	156 632
	AU version	156 633
	JP version	156 634
10	Stand for PIPETBOY (acu or comfort) and VACUBOY	155 065

Technical Data

Dimensions (LxHxW):	PIPETBOY comfort:	130 x 145 x 35 mm
	Charging station:	175 x 75 x 80 mm
	Mains adapter:	75 x 45 x 45 mm
Weight:	PIPETBOY comfort:	190 g
	Charging station:	240 g
	Mains adapter:	190 – 300 g
Batteries:	rechargeable, 2 x 1.2 V AAA NiMH/1000 mAh	
Mains adapter output:	9 VDC/500 mA	
Operating condition:	Operation:	5 – 40° C, max. 80% RH
	Storage:	-10 – 50° C, max. 95% RH

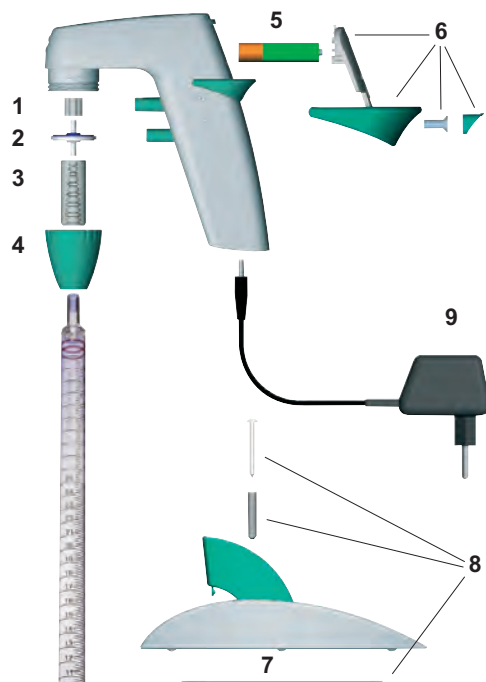

classic
156 000



volcano
156 001



pacific
156 002



royal
156 004



sahara
156 003



INTEGRA Biosciences AG
 CH-7205 Zizers, Switzerland
 Phone: +41 81 286 95 30
 Fax: +41 81 286 95 33
 info@integra-biosciences.com
 www.integra-biosciences.com



INTEGRA