

INTEGRA



MINI 96

96 Channel Portable Electronic Pipette
Affordable and Lightweight



Full or partial plates

4 models
for 0.5 to 1250 μ l pipetting

96 to 384
well microplates

AFFORDABLE AND LIGHTWEIGHT



Carry handles make it easy to move the lightweight instrument anywhere in the lab – including inside laminar flow cabinets.

A large display and comfortable grip make the MINI 96 easy to operate – no training required at all.

Full or partial 96 well plates, as well as 384 well plates, can be filled quickly and easily.

MINI 96

96 channel Portable Electronic Pipette

The **MINI 96** is the perfect solution for laboratories that want to fill 96 or 384 well plates faster and more precisely, while keeping under budget. As the most affordable 96 channel option on the market, it enables simultaneous whole plate pipetting – as well as partial plate filling – to ensure that your lab's productivity skyrockets compared to using 8 and 12 channel pipettes. Its very small size and portability allow the instrument to be placed anywhere in the lab, and it's extremely easy to use – simply turn it on and start pipetting!

Access to more productivity for every lab, every work space, and every budget.

Did you know?

The **MINI 96** is the perfect size for use in a laminar flow cabinet – and can be easily moved in and out as required.



Versatility

The **MINI 96** is available in four volume ranges to suit your application needs. The pipette's small size makes it easy to have two instruments side-by-side to perform different steps in the same workflow.



0.5 – 12.5 µl
MINI 96 Pipette


5 – 125 µl
MINI 96 Pipette

10 – 300 µl
MINI 96 Pipette

50 – 1250 µl
MINI 96 Pipette

Affordability

The **MINI 96** is simply the most affordable 96 channel electronic pipette available today. This puts significantly increased pipetting productivity within the reach of virtually every lab, freeing up time and resources for other activities.

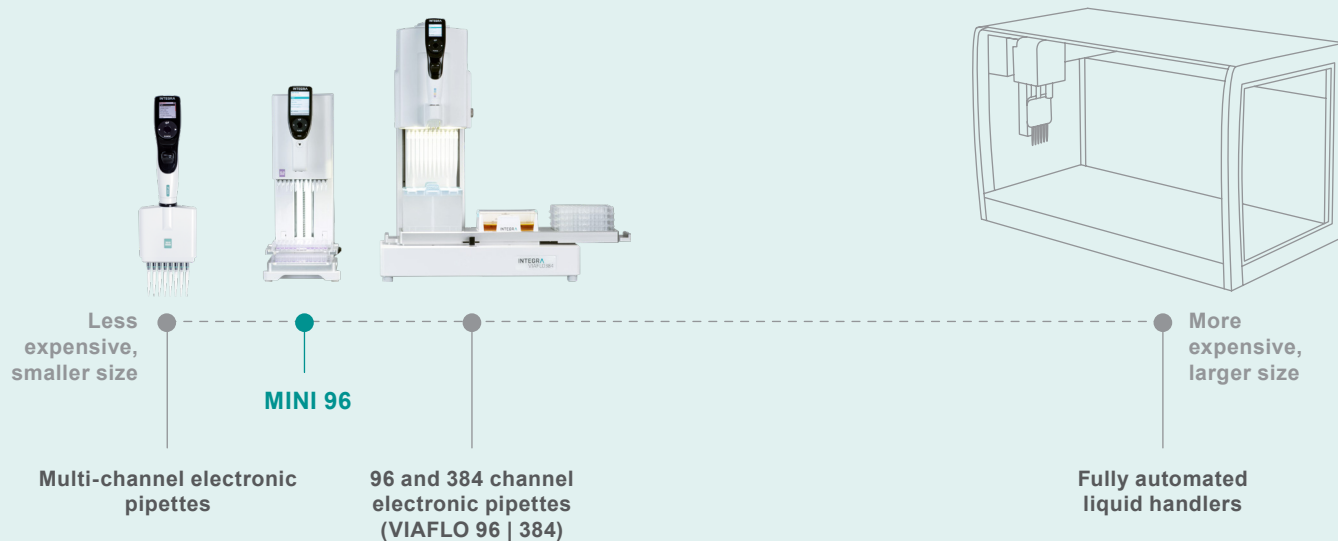

Productivity

The **MINI 96** provides a big increase in productivity – significantly faster than laboriously filling microplates with an 8 channel pipette – since an entire plate can be filled at once. This also helps to improve reproducibility and accuracy, as all channels are pipetted at the same time and speed.

Reservoir to plate transfers

Repeat Dispense mode allows ultra-fast and precise filling of multiple microplates with the same liquid, while Mixing mode ensures that dispensed liquids are homogeneous – with no strain on your hand.



Reservoir to 96 well plate transfers




Reservoir to 384 well plate transfers

Plate to plate transfers

Quick and simple duplication or reformatting of plates, plus easy supernatant removal. Aspiration heights can be precisely defined to ensure optimal pipetting results while avoiding the risk of tip crashing.



96 well plate to 96 well plate transfers

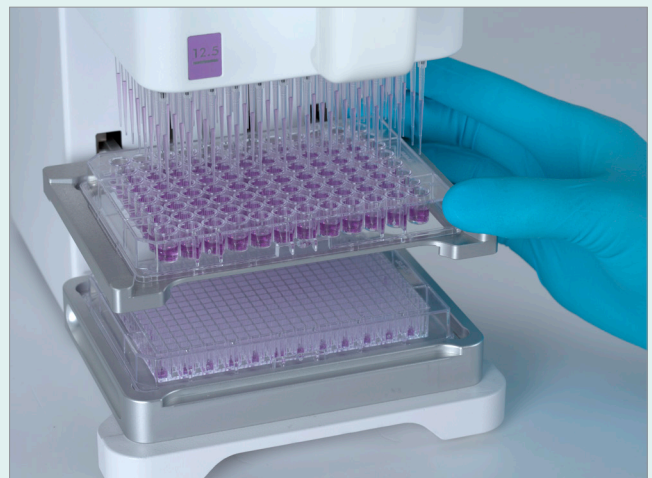

96 well plate to 384 well plate transfers

Fill half plates or single columns

Filling partial plates or performing serial dilutions is easy using the optional two position stage.



Flexibility to fill as many columns as needed...



... or perform sample transfers and dilutions within a plate.

Add a second stage



A removable second stage comes standard with each **MINI 96** to reduce labware handling. The source plate or reservoir can be placed on the base stage, then target plates can be loaded onto the second stage for filling.

GRIP TIPS®

MINI 96 features INTEGRA's advanced GRIP TIPS system. GRIP TIPS are guaranteed never to lead, loosen or fall off during a protocol, and will always attach at the exact same height, ensuring error-free pipetting. **MINI 96** even offers power-assisted tip loading and ejection to minimize the risk of fatigue or repetitive strain injuries.



💡 Did you know?




Streamline magnetic bead purifications

MAG and HEATMAG modules are user-friendly solutions for efficient magnetic bead purification. In combination with the **MINI 96**, productivity can be boosted by cleaning-up 96 samples simultaneously. The modules strong magnets move vertically by pushing a button. This enables fast bead collection with minimal bead loss for high yields.

Portability

The **MINI 96** has an ultra-compact footprint and lightweight design, making it easy to move around the lab. This makes it well suited to busy labs where multiple users are working in different places, as they can simply pick up the instrument and move it to where it's needed most – even inside laminar flow cabinets.

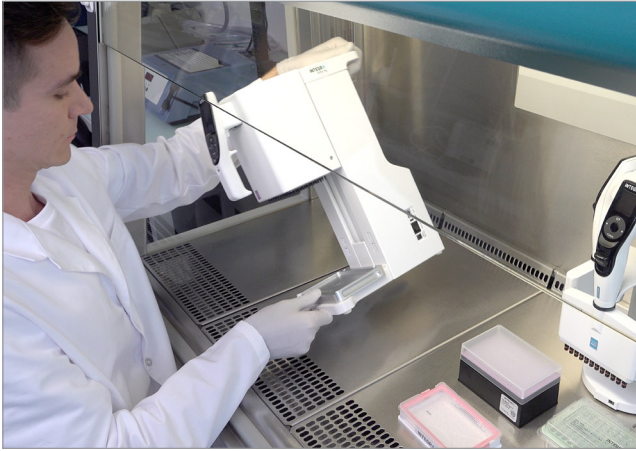

Easily moved into and out of laminar flow cabinets



Ease of use

MINI 96 is as straightforward to use as a single channel electronic pipette. It does not require any special training to operate – with simple on-screen tutorials for new users – allowing anyone in the lab to simply walk up and start pipetting within minutes of turning the instrument on. The instrument's large touch wheel-controlled graphical user interface intuitively guides the user through a choice of predefined pipetting tasks, or you can develop your own custom workflows. Together with motor-assisted operation for precise tip positioning, this makes the **MINI 96** the ideal solution for repetitive liquid handling tasks.



Ordering information

Visit the *parts and numbers* section on our website to view all available accessories and consumables.

Description	Part no.
MINI 96 Pipette, 0.5 – 12.5 µl	4801
MINI 96 Pipette, 5 – 125 µl	4802
MINI 96 Pipette, 10 – 300 µl	4803
MINI 96 Pipette, 50 – 1250 µl	4804

Description	Part no.
Two Position Stage required for loading less than 96 tips	4850
MAG module for magnetic separation	4900
HEATMAG module for heating and magnetic separation	4901
Adapter for 96 well PCR plate (MAG/ HEATMAG)	4906

GRIP TIPS for use with MINI 96

Description	Part no.
For part no. please refer to our website (www.integra-biosciences.com/GRIP TIPS)	



Automation Friendly Reagent Reservoirs

Description	Part no.	
Reservoir bases	150 ml, 8 reservoir bases, non-sterile, automation friendly (microplate format) 300 ml, 8 reservoir bases, non-sterile, automation friendly (microplate format)	6301 6305
Reservoir inserts	150 ml, 30 individually sealed, sterile, incl. 1 base and 1 lid 150 ml, 100 bulk, sterile, incl. 1 base and 1 lid 300 ml, 30 individually sealed, sterile, incl. 1 base and 1 lid 300 ml, 100 bulk, sterile, incl. 1 base and 1 lid	6317 6318 6327 6328
Accessories	PCR 96 well Cooling Block PCR 384 well Cooling Block	6250 6255


Technical data

Part no.	4801	4802	4803	4804
Volume range	0.5 – 12.5 µl	5 – 125 µl	10 – 300 µl	50 – 1250 µl
Dimensions (W x D x H)	16 cm x 26 cm x 40 cm (6.3" x 10.2" x 15.8")		16 cm x 26 cm x 44 cm (6.3" x 10.2" x 17.3")	
Weight	8.8 kg (19.5 lbs)		9.4 kg (20.7 lbs)	
Channels	96			
Compatible plate formats	96 and 384 well microplates, shallow and deep well			
Pipetting technology	Air displacement			
Operating conditions	Temperature 5 °C – 40 °C, Humidity: < 85 % non-condensing			
Power requirements	100 – 240 VAC, 50/60 Hz			

Pipetting specifications

Accuracy	± 1.5 %	± 1.5 %	± 1.2 %	± 1.8 %
Precision	≤ 1.0 %	≤ 1.0 %	≤ 0.8 %	≤ 0.8 %

Specifications presented for nominal volume. Please refer to our website for a complete list of specifications.



Contact us:

