


Reagent reservoirs: Automation friendly reagent reservoirs

Laboratory equipment | Labware | Reagent Reservoirs

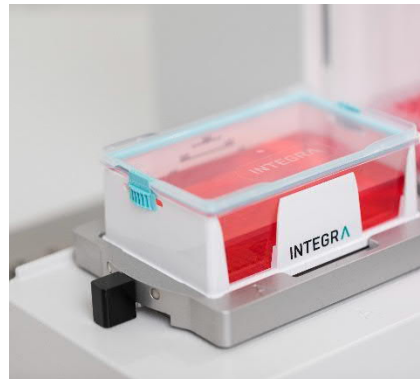

Automation Friendly Reagent Reservoirs Offer Maximum Liquid Recovery

Product description:

INTEGRA's automation friendly reagent reservoirs (150 ml and 300 ml) with standard SBS-footprint are fully compatible with all liquid handling systems including the VIAFLO 96 and VIAFLO 384 Handheld Multichannel Electronic Pipette.

The clear polystyrene reservoirs are both economically and environmentally friendly because users can reuse the reservoir base and dispose of the reservoir inserts. These reagent reservoirs are also designed to nest within one another without sticking and come in compact, space-saving packages.

In addition, the reusable base feature integrated volume graduations for accurate measurements and dual viewing windows to safely position the tips at the reservoir bottom. The anti-sealing pattern of the flat bottom, which prevents pipette tips from sealing off, and a specially formulated surface treatment result in a dead volume of less than 3 ml.

Features and benefits:

- Dual viewing windows to safely position tips directly at the bottom.
- Low dead volume due to the SureFlo™ anti-sealing array and the surface treatment of the reservoirs.
- Latching lid with gasket to prevent spillage or evaporation during short term storage.
- Less plastic waste, clear polystyrene reservoir inserts nest into a reusable base.