


## Reagent reservoirs: Multichannel pipette reagent reservoirs

Laboratory equipment | Labware | Reagent Reservoirs

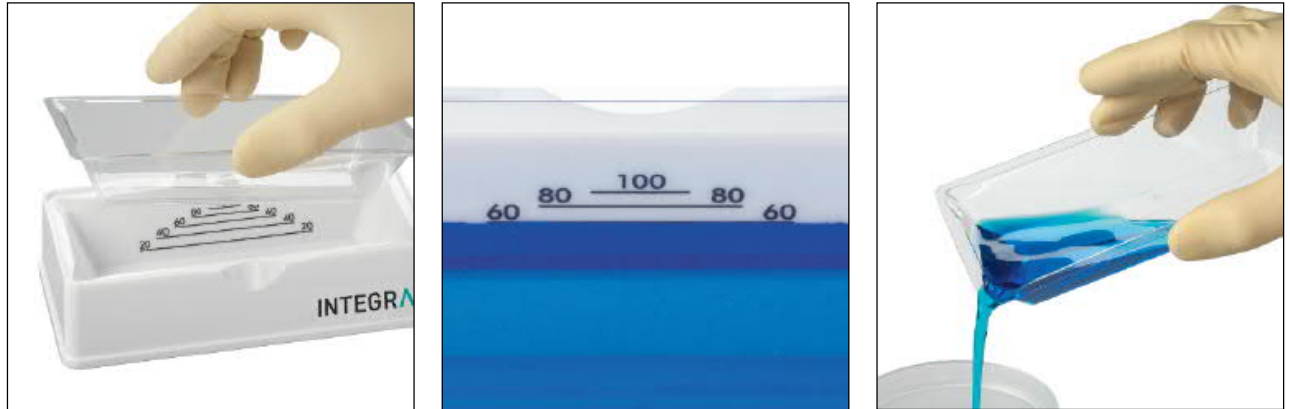

### Disposable reagent reservoirs with clearly visible volume markings

#### Product description:

INTEGRA reagent reservoirs offer a convenient solution for the temporary storage of non-solvent based liquids during any pipetting application. The system consists of disposable inserts that fit into reusable bases with clearly visible volume markings. Inserts are available as crystal clear polystyrene reservoirs or polypropylene reservoirs for improved chemical compatibility. The pour-back spouts in the corners of the reservoir allow pouring back the liquid into the original container without spills. The reagent reservoirs in sizes 10, 25 or 100 ml are stackable for space saving storage. The reservoirs are also available with the SureFlo™ anti-sealing array and hydrophilic surface treatment providing especially low dead volumes.

#### Features and benefits:

- Light is refracted so that graduations below the liquid surface become invisible, making it easier to identify the exact desired volume line.
- Pour-back spouts molded into the corners of each reagent reservoir control liquid flow and enable the easy, spill-free return of reagent to its original container.
- Base accommodates two reservoirs, one acting as a lid.
- Low dead volumes due to SureFlo™ anti-sealing pattern, inert hydrophilic surface treatment.